BEAUTIFLY

Body

Body

series

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Anti-cellulite
massager

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User manual

[beautifly.eu](https://www.beautifly.eu)

[@beautifly.eu](https://www.instagram.com/beautifly.eu)

[beautifly.eu](https://www.beautifly.eu)


[/beautiflyeu](https://www.youtube.com/beautiflyeu)

[@beautifly.eu](https://www.tiktok.com/@beautifly.eu)

f1

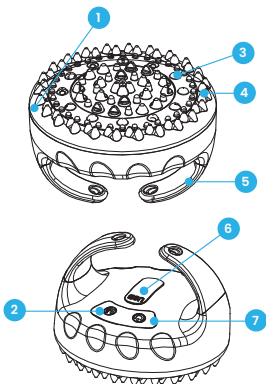f2

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- www.beautifly.eu
- www.youtube.com/beautiflyeu

Working Voltage:	3.7V
Product Dimension:	110*110x78.8 (L*W*H)
Product weight:	~198g
Material:	ABS, medical silicone
Landing:	USB

Security:

Before using the product, read this user manual and follow the recommendations contained therein when using the device.

Warning:

Do not use while charging. Do not disassemble or reassemble the product. After use, put in a safe place and keep away from children.

Device description [Page 2-f1]:

1. Massage zone
2. Mode switch
3. Massage wheel
4. Silicone massage protrusions
5. Handle
6. USB input cover
7. Power switch

Landing:

Open the USB cover and connect the cable. Connecting the cable with plug, it will start charging.

Functions:

1. Firming massage
2. Skin stimulation using special silicone protrusions.
3. Care and enhanced slimming effect thanks to photon lights.

How to use:

Waist: Massage upwards.

Belly: Massage in circular motions.

Hips and buttocks: massage in circular movements, upwards.

Thighs: massage in circular movements, upwards.

Shoulder: Massage upwards.

Contraindications:

1. Not suitable for pregnant women's skin or inflamed skin.
2. Keep out of reach of children.
3. Do not use silicone or granular emulsion with the product, so as not to destroy the massage surface.
4. If you have any bad skin symptoms, contact your doctor.

How the Body device works:

1. The device has 3 programs:

- for red light
- for green light
- for blue light

2. The device cannot be turned on while charging. The power indicator is flashing to indicate charging. When the power indicator stays on, charging is complete.

3. The device will automatically turn off after 15 minutes in working mode. Once turned on again, it can be used continuously. Continuously flashing power indicator indicates low power required for charging.

4. The device has a memory function. After turning on, the device remains in the last used mode.

5. Turn on by long pressing the switch for 2 seconds. When working in M mode, 3 modes can be switched cyclically.

Features of font lights:

1. Red light: reaches fibroblasts, which are responsible for the synthesis of collagen, elastin and hyaluronic acid. It enhances cell activity, accelerates metabolism, stimulates the skin to produce collagen proteins and delays the death of existing collagen. It also has a detoxifying and draining effect and increases the flow of the lymphatic system.

2. Blue light: inhibits the development of inflammation. Accelerates the healing of wounds, bruises and bruises. They are used against actinic keratosis.

3. Green light: works at the level of the epidermis. Indicated for reducing discoloration and sun spots. Penetrates the basic layers of the epidermis, stimulates nutrients, and helps regenerate mature skin.

More information about the company and products [page 2 - f2]:

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled.

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/UE of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/UE of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/UE of the European Parliament and of the Council as regards the list of restricted substances.

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