BEAUTIFLY

Glossy Acne

Young
series

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
f1

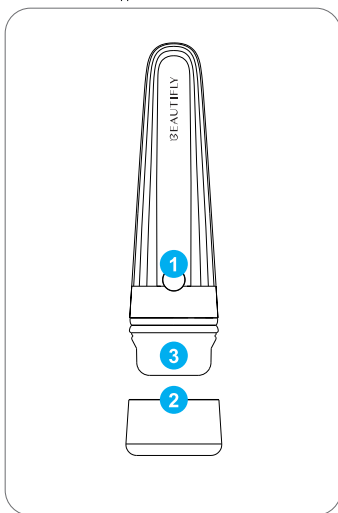f2

Size:	161x161x44.5mm
Libra:	110g
Rated power:	5W
Rated current:	1A
Battery:	3.7V polymer lithium ion
Battery Capacity:	500mAh
Included:	satin bag, USB cable, user manual

Security:

Before using the product, read this user manual and follow the recommendations contained therein when using the device.

Device description [page 2-f1]:

1. Power button
2. Protective cover
3. Treatment head

Warnings and precautions:

- Do not place or store the device near a water source such as a sink or bathtub.
- Do not immerse the device under water.
- Do not leave the device connected to the network unattended.
- Do not use the device if the cord or power supply is damaged or not functioning properly.
- Use the device only with the included power cord and power supply.
- Keep the device out of the reach of children.
- Make sure your hands are completely dry when operating the device.
- Turn off the device after use or before cleaning and maintenance.
- According to the instructions, use appropriate care creams during the treatment.

Safety rules for use:

1. Use the device only for its intended purpose as described in this manual.
2. Do not use the device on the eyelids, closed eyes, mouth, ears, thyroid gland area, breasts or genitals. Do not use the device around the eyes or mouth.
3. If you feel any discomfort while using the device, stop using it immediately.

How to use:

1. Before first use, make sure the battery is fully charged.
2. When charging, please insert the electric wire into the charger plug and then connect it to the general household power supply. To charge the product, insert the other end of the electric cable into the charging hole in the device. While charging, the LED indicator emits light and turns off when the battery is fully charged.

The charging time of the device is approximately 1.5 hours and the operating time is approximately 2.5 hours.

Power Button - On/Off: Press the POWER button to turn on the device and activate blue photon light mode.

Press the POWER button to turn off the device

Low battery information: When the device's battery is low, the red light next to the button flashes every 2 seconds.

Timing function: The timing function is 3 minutes (the signal appears once every 5 seconds).

How to use:

1. Before the treatment, clean your skin thoroughly.
2. Press the power button to turn on the device.
3. Point the device at the acne-affected area and scan it for five seconds. Move the device over the entire treatment area, making sure all acne-affected parts are scanned in turn.

4. Press the power button to turn off the device.
5. Take an appropriate amount of skin care cream and apply it evenly on the acne-affected areas. Gently massage the product all over your face.
6. Clean the treatment head of the device with a damp cloth or paper towel and apply the protective cap.

Note:

It is recommended to perform a patch test 24 hours before first use to ensure that no adverse reactions occur. Do not use the device on or around the eyes.

Recommendations for using the product:

- The treatment for the entire face lasts 3 minutes. To shorten the procedure, focus on specific problem areas. Treatment time is different for each person.
- To achieve optimal results, use the device 1-2 times daily for 1-2 weeks.
- To achieve long-lasting results after completing the full series of treatments, it is recommended to use the device 2-4 times a month or as needed.

Warnings:

- Do not use the device if you have an implanted pacemaker or internal defibrillator.
- Do not use the device if you have metal implants such as gold threads in the treatment area. However, this does not apply to fillings or dental implants.
- Do not use the device if you have been diagnosed with cancer, have had cancer in the past, especially skin cancer, or have precancerous lesions.
- Do not use the device if you suffer from serious medical conditions, such as heart problems.
- The device should not be used by pregnant or breastfeeding women.
- Do not use the device if you suffer from immunosuppressive diseases, such as AIDS, or are taking immunosuppressive medications.
- Do not use the device if you have any heat-stimulated disease in the treatment area, such as recurrent herpes simplex.
- Do not use the device if you have any active skin conditions within the treatment area such as (but not limited to) ulcers, psoriasis, eczema or rash.
- Do not use the device if you have a history of keloid scarring problem of improper wound healing or sensitive skin.
- Do not use the device if you have had surgery, laser resurfacing or deep chemical peeling in the treatment area within the last three months, or are still undergoing treatment.
- Do not use the device if you have permanent fillers in the treatment area.
- If you have temporary dermal fillers in the treatment area, consult your doctor before using the device.
- Do not use the device if Botox has been administered to the treatment area within the last three months.
- Do not use the device if you have any suspicious lesions in the treatment area.
- Do not use the device if you are taking medications, herbal remedies, vitamins or dietary supplements that may cause your skin to become dry or tender.

More information about the company and products [page 2 - f2]:

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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