BEAUTIFLY

Scrub Aqua Glow

PRO
series

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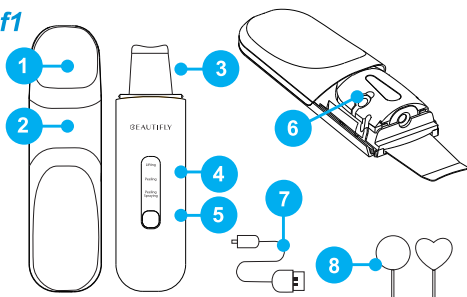
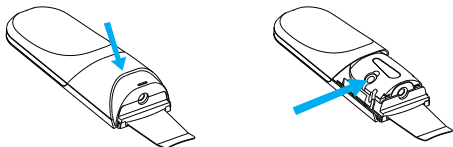
f4

f5

Ultrasound frequency:	25 kHz
Power:	2W
Output voltage:	3,7V
Power supply/charging:	battery / USB wire
Full charge time:	3,5h
Moisture protection:	IPX5

Safety

Read this manual before using the product and follow the instructions in it when using the device.

Device components [see page 2, f1]

1. Applicator lid
2. Water tank lid
3. Applicator
4. Mode and charging indicators
5. On/off button and mode control
6. Water inlet with a plug
7. USB cable
8. Pipettes for water

Modes

1. SPRAYING mode — water spraying. This mode acts as a facial moisturizer. The device turns the water into a light haze and does not leave any drops on the face. The water breaks down into small particles that are smaller than the pores on the skin, such particles quickly penetrate into the deeper layers of the skin, moisturizing it.

2. PEELING + SPRAYING mode — ultrasonic peeling + water spraying Mode for ultrasonic peeling. In this mode, ultrasound and the water spray function work simultaneously. Ultrasonic peeling is an effective stimulating facial treatment for skin care at any age. The skin is being gently and slowly exfoliated and cleansed without stretching or mechanical damaging. In just a few sessions acne disappears, sebaceous glands function becomes normalized, the upper layer of the epidermis is renewed and the skin becomes noticeably softer and more elastic. The water spray function together with the peeling keeps the skin moisturized during the procedure.

3. PEELING mode — ultrasonic peeling Mode for ultrasonic peeling without water spray function.

4. LIFTING mode — sonophoresis/micromassage. Ultrasound makes the preparations more active. It improves the process of penetration of valuable nutrients into the deepest layers of the skin, thereby multiplying the effectiveness of cosmetic products. An ultrasound massage tones muscles, massages tissues at the cellular level, visibly rejuvenates the skin, and tightens facial contours.

Device charging:

- Make sure to charge the device fully before the first use. Full charging of the device usually takes 3,5 hours (when it is completely discharged).
- Make sure that the device is dry.
- To charge the device, use the usb cable supplied with the device.
- While charging, the first indicator will blink, then the second, then the third. When the device is fully charged, indicators stop blinking.
- When the device discharges, the indicator light will blink and then the device will turn off.
- The duration of charging may vary depending on operating time and conditions.
- During Long-term storage, the device will be discharged even without use. In this case, it is necessary to recharge the device at Least monthly.
- The battery charge lasts for 1,5 hour of continuous operation.
- Warning! The device is not usable when being charged.

Control

- Modes are controlled with one button.
- To turn the device on and off, press and hold the button for 2 seconds.
- To select the necessary mode, switch the button, the indicator light will tell you which mode

is selected.

Fill the tank with water

To use mode 1 and/or 2, you must fill the tank as follows:

- Remove the applicator lid.
- Then remove the water tank Lid. For this, slide the tank lid with your thumb toward the applicator until it clicks and simply remove it.
- Remove the cap and fill with water or liquid intended for cavitation with a special pipette. Close the plug.
- Place the lid of the tank back on the device.

Warning! Do not fill the tank with Liquids other than water to avoid damage to the device. see page 2, f2]

*During the operation of the device an acoustic signal is emitted every 3 minutes. This means that it is recommended to move on to the treatment of the next facial area.

Preparation

- Before starting the procedure, it is necessary to remove impurities and makeup from the skin of the face and décollete.
- To increase the effectiveness of the ultrasound procedure, you can prepare the skin with a Light scrub.
- Work should be performed using a water-based gel/serum/lotion designed for ultrasound procedures. The main function of using cosmetics is that ultrasound increases the activity of drugs and enhances the penetration of valuable nutrients into the deepest layers of the skin. An important parameter is that the gel / serum / Lotion is an ideal conductor of ultrasonic waves.
- Perform an ultrasonic cleaning or sonophoresis procedure with progressive movements along the massage lines of the face, neck and décollete.
- After the procedure, it will be helpful to apply a face mask suitable for the skin type.

How to use?

1. SPRAYING mode — water spraying

- Press and hold the power button for 2 seconds. The device will turn on and start spraying water.
- The indicator Light on this mode will blink.

• Direct a stream of water vapor onto your skin for a few seconds to moisturize it.

*If the spray is weak, there is not enough water, add water if necessary.

Warning! Do not fill the tank with liquids other than water to avoid damage to the device.

2. PEELING + SPRAYING mode — ultrasonic peeling + water spraying

Mode for ultrasonic peeling. In this mode, ultrasound and the water spray function work simultaneously. The water sprays for 15 seconds, then stops for 5 seconds, then sprays again, etc.

*If the spray is weak, there is not enough water, add water if necessary.

- Press the power button to select this mode.
- The indicator light in this mode will not blink, but will be on and there will be a continuous us-treatment

• Perform facial cleansing in the direction of the massage.

*Proper positioning of the device applicator during the ultrasonic facial cleansing procedure.

The applicator is at an angle of 35-45° [see page 2, f3]

3. PEELING mode — ultrasonic peeling

Mode for ultrasonic peeling without water spray function.

- Press the power button to select this mode.
- The indicator Light in this mode will be on and there will be a continuous us-treatment.
- Perform facial cleansing in the direction of the massage. [see page 2, f3]

* Proper positioning of the device applicator during the ultrasonic facial cleansing procedure. The applicator is at an angle of 35-45°

4. LIFTING mode — sonophoresis/micromassage

- Press the power button to select this mode.
- The indicator Light in this mode will be on and the ultrasound will go into impulse function.
- Turn the device with the other side of the applicator and perform the sonophoresis procedure

on the face in the direction of the massage.

- Lines shown in the illustration. [see page 2, f4]

* Proper positioning of the device applicator during the sonophoresis procedure. The applicator is positioned backwards at a 45° angle.

Warning

- Do not use this device for purposes other than those mentioned in the user manual.
- Turn off the device if it is not used. Do not use the device for a long time in one place. Moving the device too fast will affect its effectiveness.
- Take away the use of the device for people with acute, chronic illness, fever, heart disease, pregnancy, facial surgery, skin rash, blood disease, during menstruation, people with pacemakers, orthodontic accommodation.
- Do not use the device around the pupil, eyelid, injured skin, mouth and intimate parts of the body.
- Please protect the device against physical impact.
- Please use the device in accordance with these instructions. Damaged device, have it repaired.
- Please keep the device away from dust, shocks, hot or excessively cold environments.
- If the packaging is damaged, before using it, please make sure that the device is working properly.
- Keep the device out of the reach of children.
- In case of malfunction, do not try to disassemble the device yourself.
- Do not submerge the device in water, do not hold it under a stream of water - it is not waterproof.
- Use only the USB cable supplied with the product.

Storage

- Keep the device dry and clean.
- Keep away from high temperature.
- Keep out of the reach of children.

Cleaning

- Rinse the peeling spatula and wipe with a damp cloth
- Do not immerse the device in water for a long time.

More information about the company and products [page 2 - f5]

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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