BEAUTIFLY

Meso Intense

Golden series

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


f1

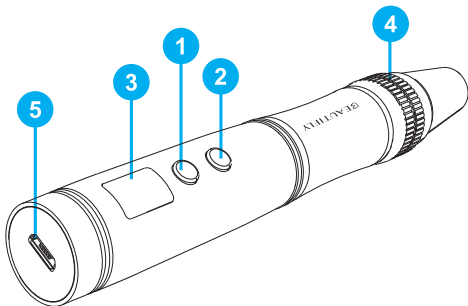f2

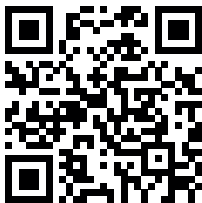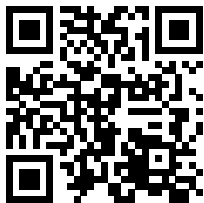Charger:	AC/DC
Entry:	AC 100-240V 50/60Hz
Exit:	DC-5V, 0.5A
Rotation speed:	10000-15000 rpm
Needle length	0-2.5mm
Adjustable speed:	6 levels
Included:	Device, charger, charging cable, user manual, 3 pcs. 12-needle cartridges

Security

Before using the product, read this manual and follow the recommendations contained therein when using the device.

Device description [see page 2, f1]

1. Power button
2. Photon light switch
3. Mode display
4. Needle length adjustment ring
5. Charging input

Installing the needle cartridge

1. Remove the needle cartridge from the sterilized packaging.
2. Align the round slot 1 with 2 of the device and push the needle cartridge in a horizontal position until it stops.
3. Turn the needle cartridge clockwise to lock it.
4. To remove the needle cartridge, turn it clockwise and push it out horizontally.

Needle depth adjustment

Adjust the depth of the needles by turning the scale adjustment ring. Turn it clockwise to shorten the needles and counterclockwise to lengthen them.

Instructions for use

1. Wash and dry your hands thoroughly before performing the treatment.
2. Use a gentle cleansing product to thoroughly cleanse your face.
3. Apply the appropriate ampoule or preparation for microneedle mesotherapy.
4. Remove the needle cartridge from the sterile packaging, attach it to the device and adjust the needle depth (the recommended needle length for first use is 0.25 mm).
5. Remove the protective cap from the needle cartridge and dispose of it in general household waste.
6. Press and hold the power button for 1.5 seconds to turn on the device.
7. Briefly press the button to switch gears. Your run information will be displayed on the screen.
8. Briefly press the phototherapy switch to select one of the 4 lights (green, red, blue, yellow).
9. Perform the treatment for 10-15 minutes, holding the device in a vertical position.
10. Gently massage in any excess cosmetics left on the skin after the treatment.

Landing

Before the first treatment, make sure the device is fully charged. Only use the USB charger included with the product.

When the power is off and while charging, the display will flash and will stop flashing

when fully charged. Do not use the charger near a TV or radio.

Cleaning and maintenance

Store the device in a vertical position to prevent cosmetics from falling into it. After use, wipe with a damp cloth.

Do not immerse the device in water or hold it under a stream of water - the device is not waterproof.

Warnings

- The needle cartridge must be moved continuously and may be held on the same area of skin for 3-5 seconds, no longer than twice.
- Temporary redness may appear on the skin after using the device. If these symptoms do not improve over time, consult your doctor.
- Avoid sun exposure, spicy foods, seafood, alcoholic drinks and • other stimulating foods for three days after treatment.
- Needle cartridges can only be used once and must be disposed of safely after use.
- Needle cartridges are for personal use only. They should not be shared with others to avoid infection.
- Do not use the needle cartridge if the disposable packaging is damaged or cracked.
- Before first use, perform a skin test on your hand. Do not use the device if your skin itches or you experience an allergic reaction.
- Parts of the device other than needle tips must not come into contact with the liquid to avoid damage to the device.
- The needle cartridge may only be used by one person.
- Do not share or reuse the same needle cartridge with others. Otherwise, it may cause skin hypersensitivity, irritation, inflammation and other problems. In the event of a breakdown, do not attempt to repair or disassemble the device yourself.

Storage

- Store the device without the needle cartridge installed.
- Store in a dry and clean place.
- Store away from high temperatures.
- Keep out of reach of children and pets.
- Place the device in a dry, ventilated place.
- Do not expose the device to direct sunlight.

More information about the company and products [page 2 - f2]

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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