

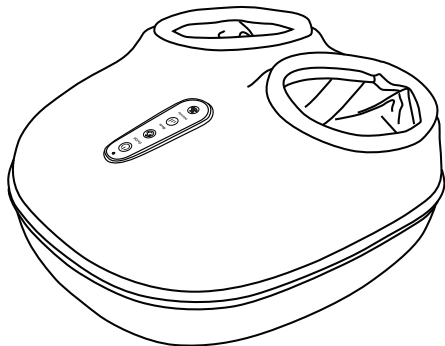
medivon.

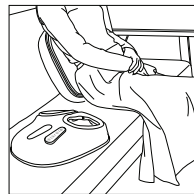User manual

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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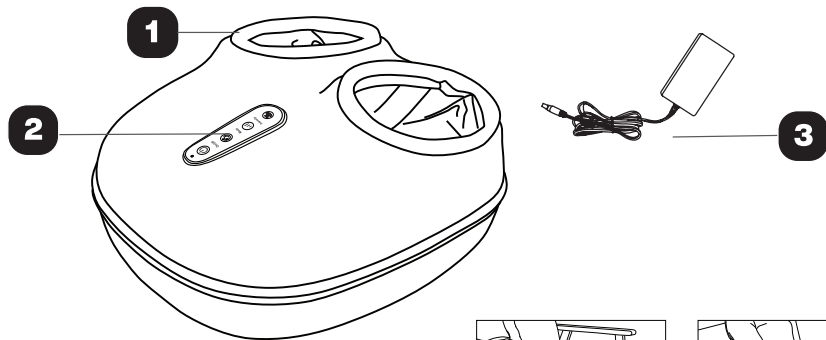1. Operation and application

1 - foot holes, 2 - remote control, 3 - adapter

The device has a rolling, kneading and heating massage function. The massager is controlled by a detachable remote control. There are three automatic massage programs with an acoustic signal. After connecting to the power supply, the device goes into standby mode.

Press the button to turn on the massager. A ready tone will appear and the device will enter mode 1 (medium massage intensity, heating on). If you want to turn off the device before the time is up, press the on/off button again; otherwise, the massager will automatically turn off after 15 minutes. After turning off the device, a long sound will appear, informing you that the massager has switched to standby.

Press the button to change one of the three available massage modes.

Press the button to change one of the three available massage intensities (low, medium, high)

Press the button to turn the heating function on or off.

There is a zipper in the center of the device. Unfold it to remove the top of the massager. The lower part of the device can then serve as a massage cushion for the back, arms and calves.

2. Dangers and warnings

People with medical implants, undergoing treatment, heart disease, high blood pressure, acute illness, pregnancy, osteoporosis, skin damage, fever above 38°C and people who are insensitive to heat should consult a doctor before use. Minors and people with intellectual disabilities may use the product only under supervision. After moving the product from a cold to a warm place, wait one hour before using it. Recommended use time is 15 minutes. Do not handle the product with wet hands and do not use it outdoors. Only use a power supply with the appropriate voltage. Do not smoke while using the product. Please use a qualified USB adapter for charging. Do not use in a humid environment. Do not use the product immediately after meals or simultaneously with other therapeutic devices. If there are any abnormalities during use, discontinue use immediately and consult a doctor. If the power cord is damaged, contact your distributor for replacement. When charging is complete or when not in use, unplug the plug from the outlet. Do not damage the power cord. Do not use on an electric heating blanket or space heater. Do not disassemble the product. In the event of a power outage, immediately turn off the product and unplug it from the power supply. Do not use the product while sleeping, intoxicated or feeling unwell. If there are any abnormalities (smoke, unpleasant odor, etc.), discontinue use immediately. Do not clean the product with gasoline, solvents or insecticides. In case of water spill, do not use the product. Do not use during a lightning storm. Before disposing of the product, please remove the battery. Keep these tips for safe and effective use of the product.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment. This device was placed on the market after 13.08.2005.

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The manufacturer hereby declares that the Medivon Artemis device complies with Directive 2014/35/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65/EU of the European Parliament and of the Council of June 8, 2011 on restrictions on the use of certain hazardous substances in electrical and electronic equipment. The full content of the declaration of conformity is available at the following internet address: www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska
Model: Medivon Artemis