BEAUTIFLY

Fresh Duo Care

Blush
series

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


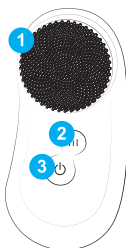
f4

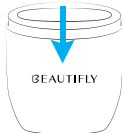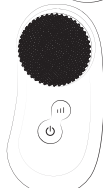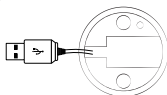f5

Device dimensions:	127.5 * 59.5 mm
Unit weight:	196,7 g
Voltage/current rating:	3.7 V 1 A
Charging time:	~ 4 hours
Water resistance:	IPX7

Safety

Please read this manual before using the product and follow the instructions in it when using the device.

Description of the unit [page 2 - f1]

1. Cleaning area / 2. Level switch / 3. Power button / mode button / 4. Smooth surface / 5. Charging station

Attention

- Keep out of reach of children.
- Do not use with clay-based facial cleansers, silicone solvents or grains to avoid damaging the soft silicone fibres.
- If you have skin or other health problems, consult a dermatologist before use.
- If you feel any discomfort during use, discontinue use immediately and consult your doctor.
- For hygienic reasons, the device should not be shared with others.
- Do not use the device or charger if they are damaged.
- Use only the charging cable supplied with the device.
- To avoid electric shock, the device should be connected to a power source equipped with a protective earthing cable.
- In skin moisturising mode, do not use products with alcohol or high oil content.
- Do not use during pregnancy, menstruation or breastfeeding.
- Do not use in case of: heart disease, severe illness, exacerbated illness, infectious disease, haemophilia, malignant tumour, allergy, psoriasis.
- Avoid using the product: on wounds, eyelids, eyeballs, plastic parts, inside the mouth, around the eyes and larynx.
- Do not rework or repair.
- Keep away from fire and heat sources.

Cleansing mode (use time: 1 minute) [page 2 - f2]

1. Remove make-up, moisten face and device. Apply facial cleanser to your face or to the cleansing area of Fresh Duo Care Blush.
 2. Hold down the Power/Mode button to activate cleaning mode. Briefly press the Level Switch to adjust the level of cleaning (you can choose from three levels of intensity).
 3. cleanse skin in circular motions.
- A low power level is recommended for first use.
4. After cleansing for 1 minute, rinse and pat your face dry. Apply skin care products appropriate to your skin condition.
- The device automatically switches off after 1 minute.

Moisturising mode (Use time: 3 minutes) [page 2 - f3]

1. cleanse your face of make-up, oil and dirt to keep your skin fresh.
 2. Before use, adjust the temperature of the tip to your skin by placing it on the back of your hand to avoid the discomfort of high temperatures.
 3. After applying skin care products to your face, press and hold the Power/Mode button to activate the skin rejuvenation mode. Briefly press the level switch button to select the intensity according to your skin tolerance (one of two intensity levels can be selected).
- A low intensity level is recommended for first time use.

4. gently glide the device over the face according to the operating technique to help the skin absorb the cosmetic.

- Absorption mode for daily skin care products: first level nourishing mode.
- Absorption mode for face masks: second level nourishing mode (sheet mask recommended).
- Do not use products containing alcohol or too much oil.
- The device will automatically switch off after 3 minutes in import mode.

Do not leave the device on one area of the face for a long time to avoid over-cleaning or skin burns.

Cleaning

1. cleaning surface of the silicone brush: clean it thoroughly after each use.

- Rinse with water.
- Dry with a tissue or dust-free towel or leave it to dry naturally.

2. smooth surface : wipe it with a damp cloth after each use.

- Note: Do not use cleaners containing alcohol, alkaline cleaners, volatile oils, solvents, nail polish remover, etc. to avoid subsequent damage to the skin.

Charging [page 2 - f4]

Battery charging time is 4 hours with a fully depleted battery, but actual charging time depends on the amount of battery charge remaining.

1. Pull the charging cable from the bottom of the charging box to the desired length and stop.

2. The unit is placed vertically in the charging box in the correct way (the corresponding function marking is indicated on the charging box).

3. Connect the power supply and place it on a flat table.

Display during charging (level lamp display)

Three horizontal lamps flash sequentially during charging.

- The first lamp flashes when the charge is not yet 1/3.
 - The first lamp is permanently on and the second lamp flashes when the charge is not yet 2/3.
 - The first and second lamps are permanently on and the third lamp flashes when the charge is not yet 3/3.
 - When all three lights are constantly on, the battery is fully charged.
 - When the battery is low, the gear indicator flashes 10 times in a row and then switches off automatically.
- 4 When charging is complete, pull out the charging plug fully and it will automatically retract.

More information about the company and products [page 2 - f5]

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled.

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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