
medivon.

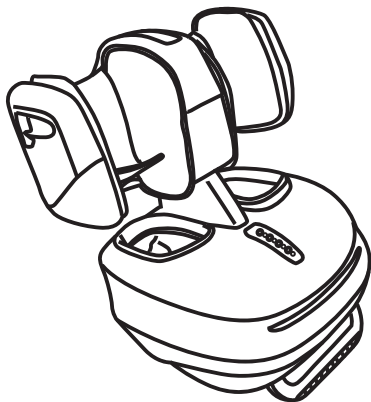User manual

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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Device Manual

After connecting the power adapter included in the set, the device is ready to work. Place your feet in the massager and press the power button to start the massager.

The upper and lower parts of the massager can be used independently. The calf/hand massager has a built-in battery that recharges while connected to the lower foot massager which is connected to the power supply.

To disconnect the parts, hold down the 'DETACH' button visible on the side of the device. Use the button only in standby mode, do not disconnect it when the device is working.

The 'MODE' button cycles through the massage modes of the device. The modes are indicated by a colored diode.

The 'STRENGTH' button cycles through the massage strength in the currently selected mode. The strength of the massage is indicated by a colored diode.

The 'HEAT' button turns the heating mode on or off

The A1 and A2 buttons are automatic modes - in case of disconnection of the parts, the automatic mode works only in the lower part of the massager.

The device must not be used by persons with injuries or who have undergone surgery in the last 3 months. Do not use the device in high temperatures and damp rooms. Absolutely keep away from water. Also, people with hypertension, hypotension or severe strokes cannot use the device. People with skin problems should consult a doctor before use. If you feel any pain or discomfort during use, stop using immediately and consult your doctor. Otherwise, it may be harmful to your health.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment. The device was placed on the market after 13.08.2005.

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The manufacturer hereby declares that the Medivon Perseus device complies with Directive 2014/30 / EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65 / EU of the European Parliament and of the Council of 8 June 2011 on restriction of the use of certain hazardous substances in electrical and electronic equipment. The full text of the declaration of conformity is available at the internet address: www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska
Model: Medivon Perseus