medivon®

User manual

MASSAGE MAT

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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1. Setup and operation

1. Secure the cushion in place by wrapping the back elastic straps over the car seat or chair.
2. Connect the adapter cable to the corresponding cable at the side of the massager.

3. Plug the home adapter into an outlet or the car adapter into your car's power receptacle.
4. Turn on the device using the controller.
5. When finished place the controller in the pouch at the side of the seat massager and unplug the adapter from the outlet or receptacle.

- Do not setup or use the massager in a bathroom or similar wet/damp areas as such an environment may reduce the lifespan of the massager mechanism.
- Ensure the massager is not exposed to direct sunlight or placed near sources of heat.


Note: This product produces a warming heat. It is not intended to be used as a therapeutic heating pad. Please allow a few minutes to reach maximum temperature. The cushion reaches a maximum temperature of 51°C (124°F) after approx. 7 minutes.

2. Safety Instructions

Please read the following instructions carefully before putting your massage device into use to ensure trouble-free operation and optimum efficiency. Please retain these operating instructions for further use!

- Do not wet, do not use pins, never remove cover.
- This item is NOT A TOY. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.
- This appliance should never be left unattended when plugged in.
- Any possible repairs may only be carried out by authorized expert staff.
- Improper use and unauthorized repairs are not permitted for safety reasons and lead to loss of warranty.
- Never touch the power plug with wet hands.
- Please avoid contact of the device with water, high temperatures and direct sunlight.
- Do not use any damaged cables, plugs or loose sockets.
- Never operate under blankets where the air opening may be blocked.
- If plugs or cords are damaged, they must be replaced by the manufacturer, a service representative or by qualified personnel.
- In case of malfunction, disconnect immediately from the mains.
- Do not use if you have skin disorders, open wounds, or swollen or inflamed areas.
- Misuse or incorrect use excludes any liability for damage.
- Do not use while you sleep.
- In order to avoid excessively stimulating the muscles and nerves, the recommended massage time should not be exceeded 60 minutes at a time.
- Each massage – even a hand massage – must be refrained from during pregnancy or if one or more of the following complaints are present in the massage area: recent injuries, thrombotic diseases, all kinds of inflammations and swellings, and cancer. It is recommended to consult a doctor before having a massage for the treatment of ailments and illnesses.
- If you are dependent on electrical aids e.g. pacemakers, please be sure to consult your doctor for medical advice before taking a massage.
- The packaging material supplied should not be used as a toy.

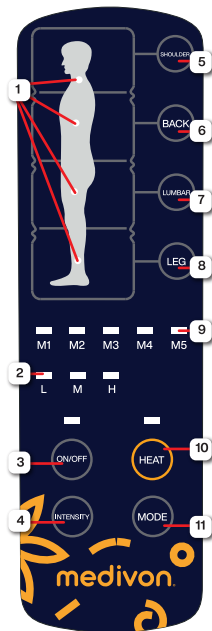3. Controller instructions

1. Body zone; indicator lights;

2. Intensity indicator lights (LO, MED, HI);

3. On/Off;

4. Intensity selector;

5. Push for Neck/Upper back massage;

6. Push for Mid-back massage;

7. Push for Lower Back/Lumbar massage;

8. Push for Leg / Thigh massage;

9. Mode indicator lights (M1 - M5);

10. Push to turn heat on / off;

11. Mode selector - choose from 5 pre-programmed massage modes.

Massage

Each massage area may be selected separately. Select high, medium or low intensity. For use at home.

Heat

Press the "HEAT" button to turn the heat function on or off. A built-in thermostat controls the temperature within a safe range. The heat may be used with or without massage.

NOTE: This product produces a warming heat. It is not intended to be used as a therapeutic heating pad. Please allow a few minutes to reach maximum temperature. Mat reaches maximum temperature of 60°C after approx. 7 minutes.

Built-in Timer

The cushion has a built-in timer that deactivates the cushion after 15 minutes. After 15 minutes of use please wait 15 minutes and simply turn cushion on again.

4. Care and Cleaning

- Disconnect from power source before cleaning.
- Do NOT place in any liquid to clean. Clean fabric cover by wiping with a damp sponge, using a mild detergent or suds. Pat dry with clean, dry cloth.
- Do NOT use abrasive cleaners, brushes, gasoline, kerosene, glass or furniture polish, paint thinner or water hotter than 49°C to clean the unit.

5. Storage

Place the unit in its box, or in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects that might cut or puncture the fabric. Do not wrap power cord around the unit. Do not hang the unit by the power cord.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment. The device was placed on the market after August 13, 2005.

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The manufacturer declares under his sole responsibility that the Cosy Vibra device complies with Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, as well as Directive 2011 of the European Parliament and of the Council /65/EU of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to the Directive of the European Parliament and 2011/65/EU as regards the list of restricted substances. The full text of the declaration of conformity is available at the internet address: www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska
Model: Cosy Vibra