BEAUTIFLY

Prestige Tri-Action

Golden

series

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
f1

f2

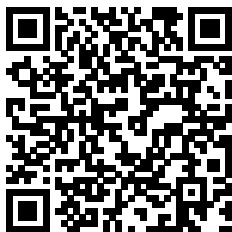Charger:	3.7V
Entry:	AC100-240V 50/60Hz
Exit:	DC5V 1A
Battery:	Built-in lithium-ion battery
Countdown timer:	15 minutes
Working time:	7-11 hours
Included:	Device, 3 attachments, charging station, power supply, user manual, warranty card

Security

Before using the product, read this manual and follow the recommendations contained therein when using the device.

Accessories [see page 2, f1]

1. Brush attachment. Suitable for scalp, face and body care. It uses electrode needles (32 pieces) and LED diodes
2. Facial attachment. Firms and smoothes facial skin. It uses electrodes and LED diodes
3. Body attachment. Firms and improves the appearance and condition of the skin. It uses an electrode section

Usage tips

- The device uses positive and negative electrodes to create a circuit that generates an electric current when it comes into contact with the skin. For best results, make sure the electrodes remain in full contact with the skin during treatment.
- The sensation caused by electrical stimulation may be different on different parts of the body. Therefore, start the treatment at level 1 (when only LED 1 is lit) and adjust the power to ensure comfort of use.
- Do not use the device unreasonably or on areas of skin with problems or abnormalities.
- After using the device, temporary redness may appear on the skin. If these symptoms do not improve over time, consult your doctor.
- Before use, make sure the device is fully charged.
- Before first use or when the battery is low, charge it using the charging station.
- Before placing the charging terminal in the charging station, carefully dry it with a cotton swab and wait until it is completely dry.
- Full charging time is approximately 2.5 hours and operating time is approximately 7 hours when using the device at level 5 (with the brush attachment, at a room temperature of 27 degrees Celsius). However, actual operating and charging times may vary depending on the previous battery level and charging conditions.

Before using the device, it is recommended to perform a patch test on the skin:

1. Apply an appropriate amount of gel or other cosmetic that conducts EMS currents to clean skin.
2. Turn on the device and test the operation of the device on the inside of the arm for approximately 1 minute, using the lowest power level (when only LED 1 is lit).
3. If you experience any unusual reactions during the test, such as rash, itching or pain, stop the procedure immediately.

Modes and functions:

The device uses two types of electrical stimulation that automatically adapt to the attachment used.

Skin lifting mode* - through impulses sent radially from the electrode diodes, it stimulates the skin and hairline muscles, providing lifting and a visible improvement in tension. Used with a brush attachment for the care of the scalp, face, arms and neckline. Used dry, however, when using a moderate level of stimulation, it can be used on damp skin or hair. It does not require prior application of gel or other cosmetics that conduct EMS currents.

Using the lifting mode on the scalp, face and body:

- Perform the treatment once a day for a maximum of 30 minutes, assuming that the session lasts from 5 to 15 minutes.
- Start with the lowest mode (level 1) to get used to the tingling sensation caused by electrical stimulation.
- Do not use the device on the nose, eyelids, ears, mouth or eyes.
- For best results, adjust the movements and level of electrostimulation to the area of the body being treated.

Head care

- Use on clean, combed hair and scalp
- Comb your hair in a slow zigzag motion, starting at the roots and working towards the back of your head
- From the back of the head, gradually move the device in a zigzag motion from above towards the opposite side of the collarbone
- Brush perpendicular to the skin, in full contact with the hairy part of the head; do not press the cap too hard to eliminate the risk of damage

Face care:

- Remove makeup and cleanse your face thoroughly
- Starting from the chin, the lower part of the jaw or the nasolabial folds, gently touch the skin with the brush, moving towards the ears; perform similar movements on the forehead, starting from its base and working towards the hairline; repeat the series 2 to 4 times

Neck, shoulder and décolleté care:

- Cleanse your skin thoroughly
- Move the brush in a smooth motion, starting from around the ear towards the neckline; from the back of the ear, go through the head towards the collarbone; finally move the device around your neck and shoulders.

*During the treatment, you will feel a tingling sensation caused by electrostimulation.

Deep EMS mode* - through a deep EMS pulse, it strengthens the muscles and increases the production of collagen and elastin. Visibly reduces fine lines and wrinkles and improves skin firmness. Stimulates blood and lymph circulation, providing the skin with a young and radiant appearance. Used with face and body pads that ensure comfort of use. It requires the use of a gel or other cosmetic that conducts EMS currents. If the skin is dry, the electrical stimulation may

be difficult to feel. It is recommended to avoid using water or lotion as they may dry out the skin and reduce the effectiveness of electrical stimulation while causing discomfort.

Using EMS mode on face and body skin:

- Perform the treatment once a day for a maximum of 30 minutes, assuming that the session lasts from 5 to 15 minutes.
- Start with the lowest mode (level 1) to get used to the tingling sensation caused by electrical stimulation.
- Do not use the device on the nose, eyelids, ears, mouth, eyes or throat.
- To ensure comfort, do not press the device too hard against your cheeks, chin or other parts of your face when using it.
- Remember that when using the device, a stronger feeling of stimulation may appear on the scalp muscles.

face care

- Before use, remove makeup and cleanse your skin thoroughly
- Use the facial care cap and apply a gel or other cosmetic that conducts EMS currents
- Starting from the chin, the lower part of the jaw or the nasolabial folds, hold the device locally for 3-5 seconds, moving towards the ears; perform similar movements on the forehead, starting from its base and working towards the hairline; repeat the series 2 to 4 times

Neck and cleavage care

- Carefully cleanse the skin and apply a gel or other cosmetic that conducts EMS currents
- Apply pointwise for 3-5 seconds, holding the device in the selected place and moving to other parts of the neck and décolleté; repeat the series 2 to 4 times

Body care:

- Carefully cleanse the skin and apply a gel or other cosmetic that conducts EMS currents
- Please note that in this mode it may be difficult to feel the EMS stimulation if your skin is dry
- Use on the muscles of your arms and legs by moving the device from bottom to top
- Apply to the belly in a zigzag motion, starting from the sides and moving towards the navel
- Use on the buttocks by moving the device from bottom to top
- Gently hold the device on the areas you particularly want to focus on
- Avoid stopping the device on the bones
- *A slight tingling sensation is felt during the treatment.

For more effective care, you can use different attachments during treatment. The lifting effect through electrical stimulation will be achieved by using a brush attachment for 15 minutes a day to improve the appearance and firmness of the skin before applying makeup. The recommended use time for each area is as follows: - Face, head, arms and neck (5 minutes).

To exercise the facial muscles and additionally improve the condition of the skin, it is recommended to repeat the treatment using the facial tip 3-4 times a week for 20 minutes before going to bed. The recommended application times for each area are as follows:

- Face, head, shoulders and neck (10-15 minutes) - Scalp (5-10 minutes) - Face, head, shoulders and neck

Please remember that the recommended time of using the device on individual parts of the body should be adapted to individual needs and the condition of the skin of the body, face and head. Please use common sense when using it and if you experience discomfort or unusual sensations, immediately stop using the device and consult a doctor.

Storage

- Keep the device dry and clean out of the reach of children.
- Keep away from high temperature.
- Store the device in a cool, dry place.

Cleaning

- For daily maintenance, use neutral detergent and a soft cloth. Do not clean the device with organic solvents or other chemicals. Do not rinse with water.
- Avoid high-pressure washing and cleaning with hard objects that may scratch its surface.
- Between treatments on different parts of the face and body, disinfect the tips using agents intended for beauty devices

Landing:

- Connect to USB for charging. It is not possible to use the device while charging.
- When the battery level is low, LED 1 quickly flashes 3 times every 3 seconds.

Precautions and safety

- Use the device in a well-ventilated room, with a temperature of 5°C-40°C and a humidity of 30%-85%.
- Keep the device away from electrical appliances, sources of high-frequency electrical radiation, etc. to avoid interference during use.
- Do not disassemble, repair or modify the device.
- When using for the first time, monitor the time (5-15 min) and pay attention to the frequency. Excessive use may cause discomfort. If you feel uncomfortable, stop using the device.
- In case of incorrect operation of the device, stop using it immediately and contact the service.

Warnings and contraindications

- If the device has not been used for a long time, wipe it down and turn it on to check whether it is working properly. Then charge your device. If any irregularities occur, contact the service.

- People with implanted electronic devices, such as pacemakers, must not use the device.
 - Using the device near a shortwave or microwave processing device (e.g. 1M) may cause the device's output to become unstable.
 - Do not use the device during pregnancy or menstruation.
 - Keep the device and its accessories out of the reach of children.
 - Do not immerse the device in any liquid and do not wipe the product with a wet cloth.
 - Turn off the device before charging.
 - Improper use or adjustment of the device may void the service warranty.
 - The device is prohibited to be used by people with wounds, fractures, scars, cancer, varicose veins and local poor skin elasticity.
 - Persons with serious heart, liver or kidney diseases or with high fever are prohibited from using the device.
- *Scalp treatment should not be performed by people with hair extensions, hair transplant recipients and people wearing certain types of hairpieces (such as adhesive tape or extensions).

More information about the company and products [page 2 - f2]

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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