BEAUTIFLY

Lumix
PRO
series

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Radio frequency instrument with
an intensively rejuvenating and
smoothing effect on the skin

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/beautifly.eu

@beautifly.eu

beautifly.eu


/beautiflyeu

@beautifly.eu

f1

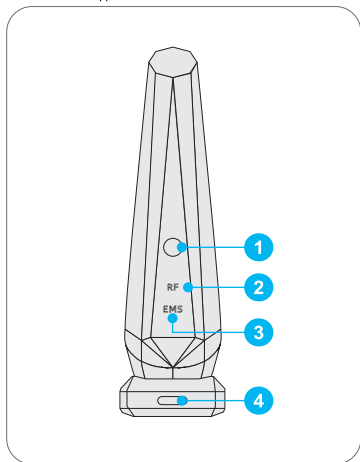f2

Size	56x187x45.8mm
Libra	Approximately 154 g
Rated power	About 10W
Rated voltage	DC 5V
Rated current	About 2A
Timepiece	6 min.
Operating mode	RF (1000KHz)/EMS 5-9μA
Battery	7.4V 600mAh

Security:

Before using the product, read this user manual and follow the recommendations contained therein when using the device.

Device description [Page 2-f1]:

1. Power button
2. RF function indicator
3. EMS function indicator
4. Charging base

Mode of action:

The Lumix PRO device, through thermal power combined with the action of EMS microcurrents, maximally stimulates collagen production for intensive rejuvenation of the facial skin.

1. Energy head plated with 24-carat gold
2. Epidermis
3. Dermis
4. The thermal power of Lumix PRO focused deep on the dermis
5. Heat passing to the deeper layers of the skin
6. Subcutaneous tissue
7. Superficial neuromuscular system

Radio waves (RF)

Radio frequency (RF) is a proven technology used in the medical and cosmetic industries. At a strictly controlled temperature, the magnetic pulse acts directly on the cells, generating thermal energy. As a result of overheating, the collagen fibers shrink, which gives an immediate effect of skin tightening. At the same time, radio frequency activates collagen production. The results are visible soon after the treatment: lifting the face contour, firming the skin, reducing wrinkles.

A series of treatments (2-3 times a week for 6-8 weeks) gives the effect of long-lasting tension of skin structures, removes wrinkles, reduces "bags" under the eyes, and restores youthful facial features.

Facial lifting massage using EMS currents (Anti-aging)

This type of massage quickly firms the face and restores skin elasticity. It smoothes wrinkles (skin ironing effect), effectively improves facial contour and reduces enlarged pores. EMS microcurrents stimulate the activity of skin cells, thanks to which the skin is intensively regenerated even at deep levels. Massage using EMS stimulates blood microcirculation, which translates into reducing discoloration, improving skin tone and maintaining proper hydration. The skin is nourished, oxygenated and full of vitality.

Energy head plated with 24-carat gold

In order to achieve the best possible non-surgical lifting effects, the head of the device is covered with 24-carat gold. Micronized gold enhances collagen stimulation, making it easier to quickly achieve lasting therapy effects. Gold in its pure form has antibacterial properties - it soothes irritations and inflammations. The energy head is a modern technological solution that makes the treatments effective and at the same time completely painless.

Skin care - reasons for performing the treatment:

1. Wash your face.
2. Apply an appropriate amount of gel/serum or other cosmetic with a slippery consistency to the area to be used (face or neck).
3. Long press the device button to turn it on.
4. Select a comfortable mode.

Important! For the first treatment session, it is recommended to set the power to low and then adjust the power level to suit your skin's tolerance for subsequent treatments.

Before using the device for the first time, it is recommended to perform a skin sensitivity test. The test should be performed 24 hours before the first use of Lumix PRO. To do this, apply a small amount of gel/serum (or other cosmetic with a slippery consistency) to the inner part of the arm and treat this area with Lumix PRO. Make sure that the skin does not react negatively to the gel/serum (cosmetics you used) or the treatment.

5. Divide your face into 4 treatment areas:

- right cheek
- left cheek
- forehead
- outer corner of the eye

6. Place the device on the skin, remember to press the device lightly to ensure full contact of the energy head with the skin.

7. Begin moving the device in very small circular motions over a small portion of the treatment area (approximately twice the size of the treatment head). The energy of the device will be released deep into the skin, providing the skin with comfort and warmth.

Important Note: Be sure to continuously move the device slowly across the skin to ensure proper, uniform energy to the treatment area. Do not hold the device in one position without moving.

8. Continue working over the entire treatment area in a slow circular motion for 2-4 minutes. Make sure to maintain uniform energy by massaging the entire treatment area.

Note: the effect of the treatment can be enhanced by performing "lifting" movements from the jawline up towards the cheeks.

9. Once you have completed the 2-4 minute treatment, you can move on to the next treatment area by repeating steps 6-8. Switch to RF + EMS mode (see "Modes" section of the manual). It is recommended to start with a low power level.

Important: Use RF + EMS mode ONLY on the jawline area. Do not use RF + EMS mode on other areas of the face.

10. Begin moving the device using very small linear movements along the jawline.

11. To ensure for optimal results, use the device for 4 minutes on each side. After the entire treatment area has been massaged, the treatment session is completed. The natural effect of the treatment may be redness of the skin (erythema) and slight swelling (swelling) in the treatment area. Both should disappear after a short time.

12. Long press the power button to turn off the device.

13. You can apply a moisturizer if necessary.

14. Clean the energy head of the device with a damp cloth or paper towel, then place the device on the protective base to protect the energy head for next use.

Warning: Do not immerse the Lumix PRO device under water.

Security mechanism:

The automatic temperature optimization mechanism in the device protects the skin against overheating.

Therapy recommendations:

- The entire face treatment lasts 8-12 minutes. To shorten the treatment time, focus only on selected areas that you want to treat. Therapy times vary for each person.
- For optimal results, use the Lumix PRO device 2-3 times a week for 6-8 weeks.
- There should be a break of at least one day between treatments on the same treatment area.
- To maintain long-lasting effects, after completing the full treatment, we recommend using Lumix PRO 2-4 times a month or as needed.

How to charge your device:

Landing

1. When using for the first time, please charge the device so that the battery is fully charged (if the device is fully charged, there is no need to charge).

2. When charging, place your device on the charging base, then insert the USB connector into the charger plug and connect it to your general household power source. The LED indicator will glow red while charging, and the red light will turn off when fully charged.

3. During the charging process, a lamp with a wavelength of 415 nm disinfects the treatment heads for 5 minutes.
4. Charging time is approximately 120 minutes and device operating time is approximately 50 minutes.

Startup and shutdown:

Press and hold the button to turn on the power in the shutdown state. The audio indicator will emit a "D" sound when the power is turned on and off.

Modes:

In boot mode, short press the button to switch successive modes from Mode 1, Mode 2, Mode 3, Mode 4, Mode 5, the buzzer will emit a "D" sound and the indicator LED color will change.

Mode 1: RF treatment (Level 1), the LED indicator will be solid green.

Mode 2: RF treatment (Level 2), the LED indicator will be solid blue.

Mode 3: EMS treatment (Level 1), the LED indicator will be solid green.

Mode 4: EMS treatment (level 2), the LED indicator will be solid blue.

Mode 5: RF and EMS treatment (LED indicator will be solid blue while automatically switching between RF and EMS therapy).

Temperature sensor- Treatment heads are equipped with a precise NTC temperature sensor.

After starting the RF treatment, the sensor will detect the skin temperature, if the temperature is lower than 42°C, the RF treatment will start working. If it is above 42°C, the RF treatment will be interrupted.

Low Battery Reminder- When the device's battery is too low, the buzzer will emit a "D" signal every 5 seconds.

Auto Timing - The auto timing function sets the time to 6 minutes. After this time, the product stops working and must be restarted to continue using it.

Safety- Please read the entire user manual carefully before using the product. Always follow all warnings and recommendations. Please observe all the following precautions when using Lumix PRO.

Warnings and recommendations:

1. Do not place or store the device near water sources such as a sink or bathtub.
2. Do not immerse the device under water.
3. After connecting the device to the power supply, do not leave it unattended.
4. Before use, check the device and its power supply for visible damage. Do not use the device if the electrical cord is damaged or if the device is dismantled.
5. Use the device only with the supplied power cord and power supply.
6. Keep the device out of the reach of children.
7. When operating the device, make sure your hands are completely dry.
8. Turn off the device after use or before cleaning and maintenance.
9. During the treatment, use a gel/serum (or other cosmetic with a slippery consistency) to ensure proper lubrication during the treatment, good RF/EMS conductivity and better skin hydration.

Safety rules:

1. The device should be used only for its intended purpose, as described in this manual.
2. Do not use the device on the eyelids, around the eyes or around the thyroid, on the lips, ears, breasts or genital area.
3. The feeling when using the device should be nice and pleasant. If you experience any discomfort, stop using the device immediately.

Do not use the device if:

1. You are under 18 years of age.
2. You have a pacemaker or internal defibrillator.
3. You have metal implants in the treatment area, such as gold threads. This does not apply to dental fillings or implants, which are not a problem.
4. You currently have or have had cancer in the past, especially skin cancer or precancerous moles.
5. You have serious comorbidities, such as cardiac disorders.
6. You are pregnant or breastfeeding.
7. You have a compromised immune system due to an immunosuppressive disease such as AIDS, or you are taking immunosuppressive drugs.
8. You have any disease simulated by heat, such as but not limited to recurrent herpes simplex (in the treated area).
9. You have a history of keloid scars, poor wound healing, or sensitive skin. You have had surgery, laser

treatment, or a deep chemical peel in the treatment area within the last three months and are still in the full healing process.

10. You have permanent fillings in the treatment area.

11. If there are temporary dermal fillers in the treatment area, consult your doctor before using the device.

12. Botox has been administered to the treatment area within the last three months.

13. You have suspicious lesions in the area being treated.

14. You are currently taking medications, herbal preparations, vitamins or dietary supplements that may increase skin sensitivity or dryness.

More information about the company and products [page 2 - f2]:

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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