BEAUTIFLY

Derma

Blush

series

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
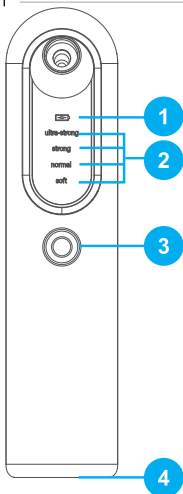
f4

f5

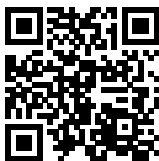Tension:	5V
Current:	450mAh
Battery:	Li-ion battery
Power:	1.6W
Size:	163.5x52x38mm
Libra:	205g
Suction power:	50-60 kPa

Safety

Read this manual before using the product and follow the instructions in it when using the device.

Part description - device [page 2 - f1]

1. Charging indicator
2. Suction power indicator
3. Power button/suction power adjustment
4. USB port

Parts description - heads [page 2 - f2]

1. Small round suction head. Low suction power For people who have dry or sensitive skin and those who start treatments with the device.
2. Oval suction head. Medium suction power/special shape. For lifting and firming treatments.
3. Round suction head. Gentle/medium suction power. For people who have oily skin, acne or blackheads.
4. Small microcrystalline suction head
5. Oval microcrystalline suction head
6. Round microcrystalline suction head

Service

On/Off

Press and hold the power button for 3 seconds, the product is turned on, the display is on, and the suction power displays 1 level (soft).

Power regulation

Short press the power button to adjust the suction power. Each short press increases the suction power by one level, in the order: soft > normal > strong > ultra strong.

Application

- Discoloration
- Scars
- Stretch marks
- Enlarged, highly visible pores
- Blackheads and acne (avoid inflammatory and purulent lesions)
- Wrinkles

- Skin lacking a healthy color with the first signs of photoaging

Contraindications

- Active purulent acne
- Rosacea
- Retinoid therapy
- Acid treatments, lasers - at least 4 weeks before microdermabrasion
- Psoriasis
- Viral, bacterial and fungal diseases of the skin
- Vascular skin
- Tendency to scarring
- Broken skin (fresh wounds)
- Skin tumors or suspicious moles in the treated area
- Fillers in the face - at least 4 weeks before performing microdermabrasion, otherwise the preparation may move.
- Do not wash the device with water, the device is not waterproof.

Microdermabrasion [page 2 - f3] - move in the direction of the arrows

Lifting massage [page 2 - f4] - move in the direction of the arrows

Set the suction head to a 90 degree position and slowly move it over your face. Do not hold in the same part for a long time. (Before removing the suction head from the skin, tilt it to the side to release the pressurized air. This will allow you to safely remove it without irritating the skin. Do not pull or move it directly to avoid pulling the skin and causing redness.)

More information about the company and products [page 2 - f5]

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled.

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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