
PL EN DE FR ES PT IT CZ SK SL EE LT LV
NB FI HU RO AR RU EL HR SV DA NL TR BG

BEAUTIFLY

BeautiLight Duo

Golden

series

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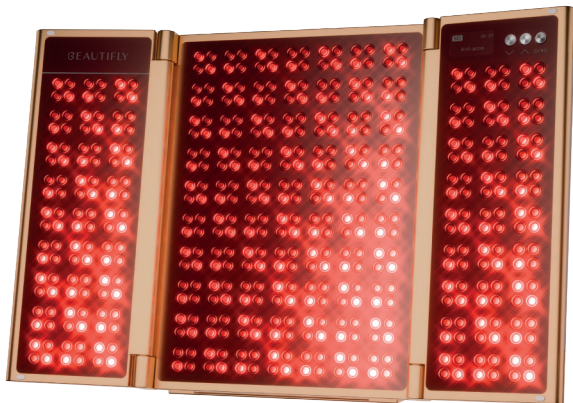Device description (see page 2, f1)

1. LED indicators
2. Display
3. Adjustment button (DOWN)
4. Adjustment button (UP)
5. Power button
6. DC power port
7. Base
8. Protective goggles
9. Power cable

Technical specification:

Dimensions: 38x27x4cm

Net weight: 2,6 kg

Power supply: AC 100-240 V / 50-60 Hz → DC 24 V / 5 A

Power consumption: up to 120W (Full Spectrum)

Wing adjustment: 0-120°

Number of light modes: 12

Number of light intensity levels: 3

Number of LEDs: 496

Wavelength spectrum: 465 nm ±10 nm (blue light), 590 nm ±10 nm (yellow light),

633 nm ±10 nm (red light), 850 nm ±10 nm (near infrared), 1064 nm ±10 nm (deep infrared)

Package contents: 1x main device, 1x protective goggles,

1x power cable, 1x user manual, 1x warranty card

Device operation

1. To turn the device on/off, open both doors or hold the power button for 1,5 s.
2. After turning on the device, you can adjust the intensity and mode using the adjustment button (DOWN) and the adjustment button (UP).
3. To change the light mode, click the adjustment button (DOWN) or the adjustment button (UP).
4. To change the light intensity level, hold the adjustment button (DOWN) or the adjustment button (UP).
5. Recommended use: Use the device 3–5 times a week, no longer than 30 minutes per day.
6. If the device is not used in standby mode, it will automatically switch off after 2.5 minutes.

Preparation for use

1. Before use, remove makeup and thoroughly cleanse the face. Do not use cosmetics containing photosensitizing or whitening ingredients.
2. Connect the device to the supplied power adapter.
3. Put on protective goggles – during irradiation the device emits intense LED light.
4. Open both wings or hold the power button for 1,5 s to turn on the device.
5. After setting the mode and energy, briefly press the power button to start the treatment. After five sound signals, irradiation begins.
6. To pause the treatment, press the power button – the time will be remembered. After 2.5 minutes of inactivity, the device will turn off automatically. The memory function saves the settings from the last use. The memory function saves the settings from the last use.

The device has 12 light modes:

1. M1 Pure Radiance Ritual (465+590+633+850+1064nm) - radiance and tone unification
- Use a gentle, non-irritating cleanser to cleanse the facial skin, apply a moisturizing skincare product without photosensitizing ingredients for basic skincare, and after absorption set the whitening mode for irradiation. After irradiation, whitening skincare products may be used. Daily attention should be paid to moisturizing and sun protection.
2. M2 Lift & Glow Therapy (633+850+1064nm) - firmness and elasticity
- Use a gentle, non-irritating cleanser to cleanse the facial skin, apply a moisturizing skincare product without photosensitizing ingredients for basic skincare, and after absorption set the firming mode for irradiation. Use regular care to achieve radiant skin.
3. M3 Clarity Light Mode (465+633nm) - balance for skin with imperfections
- Use a gentle, non-irritating cleanser to cleanse the facial skin, apply a moisturizing skincare product without photosensitizing ingredients for basic skincare, and after absorption set the anti-acne mode for irradiation. Good lifestyle habits (healthy diet and adequate sleep) are also essential for skin protection; it is recommended to go to bed early daily and avoid high-sugar products and spicy foods.
Note: Do not use in case of skin inflammation or while taking medications.
4. M4 Calm Repair Ritual (590+633nm) - soothing and regeneration
- Cleanse the facial skin with a gentle, non-irritating cleanser, use moisturizing skincare products without alcohol and photosensitizing ingredients for basic skincare, and after absorption set the regeneration mode for irradiation. The skin barrier is fragile during sensitivity periods, so daily attention should be paid to moisturizing and sun protection.
5. M5 Harmony Balance Light (465+590+633nm) - sebum control

- Use a gentle, non-irritating cleanser to cleanse the facial skin, apply a moisturizing skincare product without photosensitizing ingredients for basic skincare, and after absorption set the balancing mode for irradiation.

Since water-oil balance disorders of the skin are closely related to lifestyle habits, a healthy diet, regular physical activity, and adequate sleep are recommended.

6. M6 Intima Care Mode (465+633+850nm) - gentle intimate care

- Clean the intimate areas and then set the intimate mode for irradiation.

Note: In case of abnormal symptoms such as rash, swelling, spots or ulcers in intimate areas, go to the hospital for medical consultation.

7. M7 Soft Underarm Glow (633+850nm) - suitable for underarm care

- Clean the underarm area and set the underarm mode for irradiation. Daily washing of the underarm area is recommended. Choose products that inhibit sweat odor and absorb moisture to effectively control sweat gland secretion and odor problems.

8. M8 Derma Restore Mode (590+633nm) - suitable after microneedling, gold microneedling, needle cleansing and other aesthetic medicine treatments requiring regeneration

- It is recommended to follow the requirements of aesthetic medicine and use skincare products without photosensitizing ingredients for moisturizing. After absorption, set the special mode for irradiation and pay daily attention to moisturizing and sun protection.

9. M9 Deep Infusion Ritual (590+633+850nm) - suitable for dry skin with water deficiency and to support essence absorption

- Cleanse the facial skin with a gentle, non-irritating cleanser, apply moisturizing skincare products without photosensitizing ingredients for basic skincare, and after absorption set the penetration mode for irradiation.

10. M10 Smooth Silhouette Therapy (633+1064nm) - suitable for postpartum care of pregnant women, improves the appearance of stretch marks

- Keep the skin clean when using the pregnancy stretch mark mode. Clean the skin promptly after sweating, do not allow sweat to remain on the skin as this may cause discomfort in the treatment area. The effect of the pregnancy stretch mark mode will be better when massaging with oil during the irradiation process.

11. M11 Scalp Vital Light (633nm) - suitable for scalp health care

- Use a gentle, non-irritating cleanser to wash hair, and after hair dries set the hair care mode for irradiation. The hair care mode helps improve scalp sebum secretion, clogged hair follicles, excessive dandruff and other scalp problems.

12. M12 Total Renewal Spectrum (465+1064nm) - suitable for multifunctional skin care

- Use a gentle, non-irritating cleanser to cleanse the facial skin, apply a moisturizing skincare product without photosensitizing ingredients for basic skincare, and after absorption set the multi-care mode for irradiation. Daily attention should be paid to moisturizing and sun protection in everyday life.

Cleaning and maintenance

1. After each use, the device should be placed in its packaging for safe storage.

2. To maintain session effectiveness, clean the LED panels after each use. You can use a soft cloth or wet wipes for this purpose.

3. Regular cleaning and proper maintenance ensure optimal operation and extend the service life of the device.

1. The device contains small parts. Keep them out of reach of children to avoid accidental swallowing.

2. Do not use the device for purposes other than cosmetic.

3. Do not use the device for animal care, as this may cause injury to them.

4. Do not use the device for purposes other than those specified in the user manual, as this may lead to bodily injury or malfunction.

5. Do not allow persons with physical or mental disabilities to use the device independently without supervision, as this may lead to accidents or bodily injury.

6. There is a risk of skin irritation or inflammation if the product is used more frequently than recommended.

7. Do not use the device if it is damaged, as this may cause accidents or injury.

8. Do not use the device near flammable substances (such as alcohol, solvents, gasoline, nail polish remover or hair spray), as this may cause fire, electric shock or bodily injury.

9. Do not charge, use or place the device near fire, heating devices or in high-temperature places, as this may cause overheating, ignition or explosion.

10. Do not use the device near water sources, in bathrooms or during bathing, as this may cause electric shock, short circuit or fire.

11. Do not charge or immerse the device in water, as this may cause electric shock, short circuit or fire.

12. Avoid dropping the device and collisions with the high-power lamp and controller, as this may cause malfunction. If the device has been damaged as a result of a fall or impact, do not use it.

13. Turn off the power before cleaning or storing the device, as this may cause electric shock, short circuit or fire.

14. Do not use the power adapter or power cable if there is any damage, heating or loosening of the adapter socket, as this may cause fire, electric shock, short circuit or burns.

15. Do not apply to areas of skin where there is an increased likelihood of skin cancer, on open wounds or rashes, on swollen, reddened, irritated, infected or inflamed skin areas or skin lesions.

16. Use only the original power adapter supplied by the manufacturer. During use, ensure that the adapter plug is directly connected to the socket. Do not use extension cords, as this may cause electric shock, short circuit or fire.

17. Do not scratch, modify or excessively bend the power cable, do not pull, twist or place heavy objects on it, as this may cause fire, electric shock, short circuit or burns.

18. For at least 48 hours before using the device and for at least 48 hours after its use, do not use tanning beds or expose the skin to prolonged sunlight.

19. During use, the AC voltage must be 100–240 V. Use the original power adapter, otherwise there is a risk of fire, electric shock, etc.

20. Do not disassemble, repair or modify the device yourself, as this may cause fire, electric shock or bodily injury. If repair is required, contact after-sales

service.

21. The device is intended for home use, in a dry environment and at an air temperature of 15 to 30°C; do not use the device in places where sunlight falls directly on the skin.
22. If abnormalities occur during use, immediately disconnect the power adapter. In case of unusual smoke, noise or odor, immediately stop use.
23. Do not use the device with flammable cosmetics, as this may cause fire, electric shock or bodily injury.
24. Do not use the device simultaneously with a pacemaker, ECG scanners or other electronic devices, as this may cause instability of medical electronic devices.
25. The device is intended for home use.
26. Failure to follow the instructions may cause eye damage (even blindness) or skin damage.
27. If you are undergoing therapy with photosensitizing drugs or traditional Chinese medicine, consult a doctor before using the device.
28. The device is not intended to treat diseases. In case of rashes, itching of the skin, fungal infections, skin nodules or fibromas, medical consultation is recommended.
29. Before use, check for visible damage to the transmission mirror or housing. In case of damage, discontinue use.
30. Before first use, perform a light sensitivity test. In case of skin hypersensitivity reactions, immediately discontinue use.
31. If any abnormalities or discomfort occur during use, immediately stop using the device.
32. Use only the original power adapter included with the product.
33. Do not use or store the power adapter near damp clothing or in humid places.
34. Do not use or store the device in a place where the bulb and LED are directed directly at each other.
35. Before daily use of the device, remove makeup and cleanse the area intended for irradiation.
36. Before use, remove metal objects, glasses and contact lenses from the irradiation area.
37. When using for the first time or after a long break, ensure that the buttons and indicators function properly.
38. After use, keep the device and power adapter dry and clean.
39. Disconnect the device from the power adapter during cleaning. After cleaning, store the device in a dry and clean place.
40. After use, the skin may be slightly warm and dry. It is recommended to use moisturizing skincare products.
41. In case of severe discomfort or unusual sensations, immediately stop use and consult a doctor.
42. Select the appropriate operating mode according to individual needs.
43. The total irradiation time of the same area in one day should not exceed 30 minutes to avoid skin damage.
44. Do not use corrosive solvents (thinner, gasoline, nail polish remover, etc.) to clean the housing, as this may damage the surface of the device.

Limitations and contraindications for using the device:

1. People with skin allergies, skin diseases, atopic dermatitis, sensitive skin or a tendency to photoallergic reactions
2. People suffering from serious heart, liver or hematopoietic system diseases, nervous system diseases, cancers, facial neuropgia, hemophilia, diseases accompanied by fever and photosensitive epilepsy
3. Pregnant and breastfeeding women
4. People with capillary diseases caused by the use of anabolic steroids or with liver function disorders
5. Areas of the body where there is a break in skin continuity, such as burns, blisters, cuts, wounds, active skin diseases or infections
6. People taking medications that may cause photoallergic reactions - the device should be used only with special caution and in accordance with the doctor's recommendations
7. Occurrence during or after the treatment of abnormal symptoms at the irradiation site, such as itching, swelling, erythema, blisters or other skin reactions - in such a case, immediately stop using the device and consult medical personnel
8. Deterioration of well-being during the treatment - immediately stop the procedure and contact medical personnel
9. The device is not intended for use by children

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
[EN] Used electrical and electronic equipment, as well as batteries and accumulators, must not be mixed with other waste. This is indicated by the symbol of a crossed-out bin placed on the product packaging. Electrical and electronic equipment, as well as batteries and accumulators, may contain hazardous substances, mixtures, and components harmful to the environment and human health. Improper handling or damage to them can lead to health risks or environmental pollution. Proper handling of used electrical and electronic equipment, as well as batteries and accumulators, in households positively contributes to reuse and recovery, including recycling of the used equipment. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and served the same functions as the equipment sold. When delivering equipment intended for households, the distributor is obliged to collect used equipment from households free of charge at the delivery location, provided that the used equipment is of the same type and served the same functions as the delivered equipment. A distributor operating a retail unit with a sales area, as defined in Article 2, point 19 of the Act of March 27, 2003, on spatial planning and development, of at least 400 m² dedicated to the sale of equipment intended for households, is obliged to accept used equipment from households free of charge at this unit or in its immediate vicinity, provided that none of the external dimensions exceeds 25 cm, without the necessity of purchasing new equipment intended for households. The product was placed on the market after August 13, 2005.

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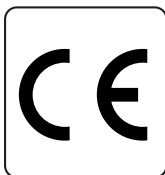[EN] Summary declaration of conformity: The manufacturer declares under its sole responsibility that the BeautiLight Duo Golden product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility; applied standards: EN IEC 55014-1:2021, EN IEC 55014-2:2021, EN IEC 61000-3-2:2019+A1:2019, EN 61000-3-3:2013+A1:2019+A2:2021; with Directive 2014/35/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits; applied standards: EN 60335-1:2012+A11:2014+A13:2017+A1:2019+A14:2019+A2:2019+A15:2021+A16:2023, EN 60335-2-27:2013+A1:2020+A2:2020, EN 62233:2008; and with Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU with regard to the list of restricted substances; applied standards: IEC 62321-3-1:2013, IEC 62321-4:2013+AMD1:2017, IEC 62321-5:2013, IEC 62321-6:2015, IEC 62321-7-1:2015, IEC 62321-7-2:2017, IEC 62321-8:2017. The full text of the Declaration of Conformity is available on the website www.beautify.eu/deklaracje

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