BEAUTIFLY

Aura Glam

PRO


series


Maska kolagenowa LED
LED Collagen Mask

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
f1

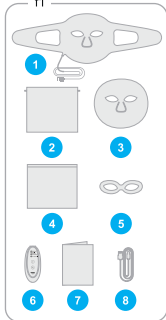f2

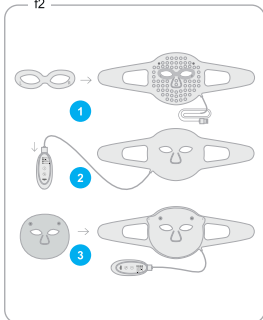f3

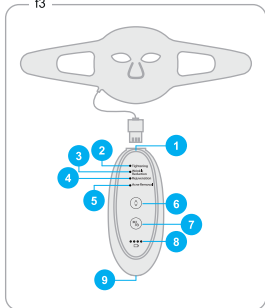f4

Device Size:	659x232x38mm
Net weight:	284g
Charging way:	USB-C charging
Charging time:	about 4 hours
Battery:	3000 mAh lithium polymer battery
Tension:	5 May

Security:

Before using the product, read this user manual and follow the recommendations contained therein when using the device.

Set contents [page 2- f1]:

Collagen LED mask (1.) Accessory bag (2.) Cooling mask (3.) Storage bag (4.) Light protection eye mask (5.) Remote control (6.) User manual (7.) Cable for USB-C charging (8.)

Configuration and first use [page 2- f2]:

1. Install the light protection glasses into the holes of the collagen mask.
2. Connect the mask with the remote control using a USB cable.
3. Instructions for using the cooling mask: After unpacking, take out the cooling mask and put it into the accessory bag. Store in the refrigerator for at least 10 minutes before using. When you need to use the cooling function, take it out of the refrigerator. To assemble the cooling mask, align the magnetic parts of the cooling mask with those placed on the collagen mask.

Remote control [page 2- f3]:

Connection to main unit (1.) Wrinkle reduction mode (2.) Firming mode (3.) Rejuvenation mode (4.) Support mode to alleviate skin imperfections (5.) Intensity change button. The remote control has three modes: low (the indicator lights up dimly), medium (the indicator flashes quickly), high (the indicator lights up constantly). (6.) On/off and mode switch button (7.) Battery indicator (8.) USB-C charging port (9.)

Light sensitivity test:

Before using the device for the first time, check whether you are allergic to light: place it on your forearm for 8 minutes. If skin becomes red or itchy, stop use and consult a dermatologist. Always clean your face before use to avoid reactions to cosmetics, and after use, don't worry if your skin becomes slightly red - this is normal and will go away after 1-2 hours. If baking, stop using the appliance immediately.

How to use:

1. Cleanse your face and dry it.
2. Fasten and adjust the headband of the LED collagen mask.
3. Press the power button to turn on the device.
4. Press the mode button to select the desired mode, then press the level button to adjust the energy level.
5. When using the cooling function, remove the cooling mask from the refrigerator and magnetically place it on the mask. After putting it on your face, activate the firming mode on the remote control.
6. After four minutes, the device will turn off automatically. Remove the mask and clean it (only clean the mask, not the electronic parts).

Charging method:

The product should not be used while charging. (Note: The device will not work and buttons will not respond while charging.) It is recommended to fully charge the device before first use. Connect the USB-C cable to the remote control, then plug the USB plug into the adapter. When the power is turned on, charging will start. When the battery indicator on the remote control starts flashing, the device is charging. When all four battery indicators stop flashing, charging is complete.

Maintenance and care of the device:

Store this product in a dry place at room temperature and avoid exposure to dust, direct sunlight or humid conditions.

Clean the mask before and after use (mask only, no electronic parts). Note: the controller and cables should not come into contact with water. After use, make sure to disconnect the remote control from the mask and from the power source.

The product contains a lithium polymer battery. When disposing of your device, please take it to your local battery recycling facility. It is recommended to separate it from other household waste to protect the environment.

Precautions:

Please keep the product away from small children to avoid treating it as a toy. Handle the product carefully when using. This device is designed as a home electronic product and is not suitable for diagnosis or treatment. It should not be used by pregnant women, preparing for pregnancy or breastfeeding, people sensitive to light or taking photosensitizing drugs, with genetically determined skin diseases or conditions such as lupus, photosensitive eczema, or albinism. Note: the controller and cables should not come into contact with water. Do not use while charging. Do not use the device if it is damaged. Do not repair, disassemble or modify the device yourself. Any irregularities should be assessed by the service. Do not drop or hit the device to avoid malfunction. After each treatment, wipe the product with a damp cloth. Do not wash the device under water or submerge it. Do not use corrosive agents to clean the device. Store it in a cool, dry place.

More information about the company and products [page 2 - f4]

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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Importer: Beautifly Sp. z o.o.

Plac Bankowy 2, 00-095 Warszawa, Polska

Model: Aura Glam