# BEAUTIFLY

Therm Body Revive

## Golden

series

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


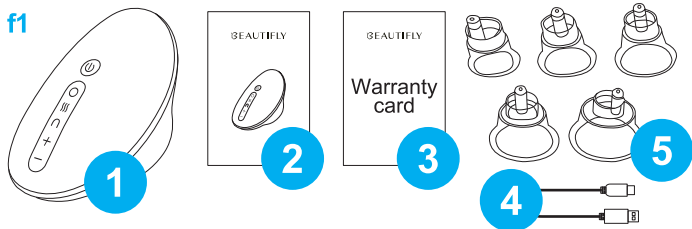
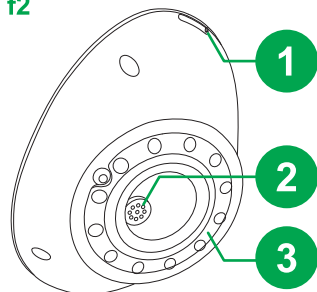
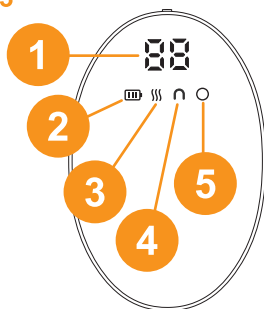
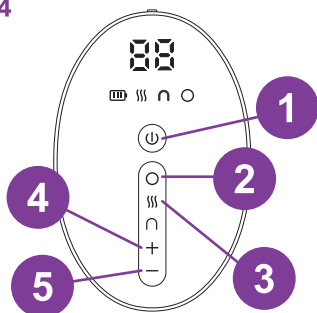
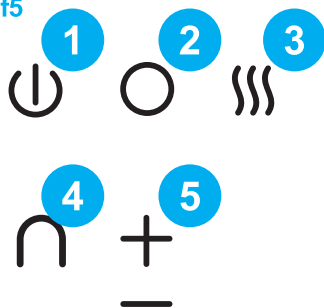
f6

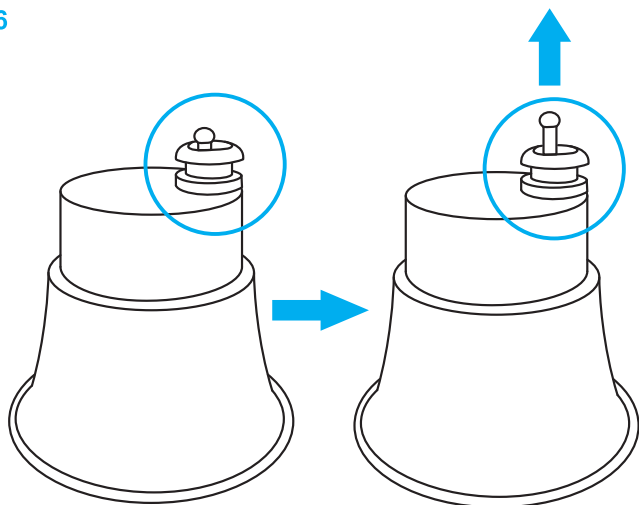f7
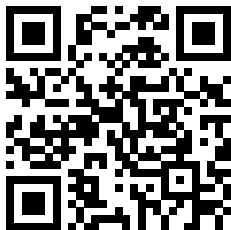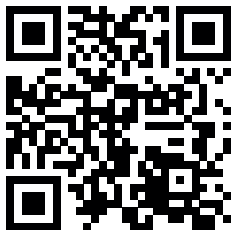
|                   |                 |
|-------------------|-----------------|
| Input voltage:    | 5V/1A           |
| Size:             | 130*86*58       |
| Libra:            | 316g            |
| Charging time:    | 4 hours         |
| Working time:     | 15 minutes      |
| Battery type:     | lithium battery |
| Battery capacity: | 2500mAh         |

#### Included [ see page 2. f1 ]

1. Device
2. User manual
3. Warranty card
4. USB charging cable
5. Vacuum massage pads

#### Device description [ see page 2. f2 ]

1. Charging port
2. Particulate filter
3. Compress

#### Panel display [see page 2. f3]

1. Numerical indicator
2. Battery indicator
3. Heating indicator
4. Pressure indicator
5. Heart rate indicator

#### Operation panel [see page 2. f4]

1. On/Off
2. Pressure button
3. Heating button
4. Intensity + button
5. Intensity button -

#### Control panel description [ see page 2. f5 ]

1. On/Off - press and hold for 3 seconds to turn on the device. A short press changes the vacuum intensity, from 1 to 5 levels.
2. Pressure Hold/Pressure Release - the first short press holds pressure, the second short press releases pressure, and the third short press creates a vacuum.
3. Heating - short press of the heating button changes the temperature from level 1 to level 3. (40°C, 45°C, 50°C)
4. Negative pressure button - press the button to turn on, press again to turn off.
5. Pressure intensity button - press the button to turn on, short press changes the intensity, between 1 and 10 levels of intensity.

#### User manual

1. Press the "on/off" button for 2 seconds to turn on the device, then short press the button to select the intensity level that suits you, from level 1 to level 5.
2. If you need to turn on the heating function, press the "Heating" button and select a level from 1 to 3 levels that is suitable for you.
3. After use, short press "Pressure/Release" and release the pressure to turn off the device.

#### Instructions - massage with accessories [see page 3. f6]

1. Press "on/off" for 2 seconds to turn on the device and select the appropriate suction cup.
2. Fit the suction cup to the device, and then select the appropriate massage intensity.
3. After use, pull back the pressure relief valve on the device and remove it.
4. Press "on/off" for 2 seconds to turn off the device.

#### Storage

- Keep the device dry and clean, out of the reach of children.
- Store away from high temperatures.
- Store the device in a cool, dry place.

#### Cleaning

- Keep the device clean. For daily maintenance, clean it with a neutral detergent and a soft cloth. Do not clean the device with organic solvents or other chemicals. Do not rinse it with water.
- Avoid high pressure washing and scratching the surface of the product with hard objects.

#### Precautions and safety

- The operating environment of the device should be well ventilated, the temperature should be kept at 5°C~40°C, and the humidity should be 30%~85%.
- Keep the device away from electrical appliances, sources of high-frequency electrical radiation, etc., to avoid interference during use.
- Under no circumstances should you disassemble, repair or modify the device.
- When using for the first time, time should be controlled (5-15 min). Moreover, pay attention to the frequency. Excessive use may cause discomfort. Stop using the device if you feel uncomfortable.
- In case of incorrect operation of the device, stop using it immediately and contact the service.

#### Warnings and contraindications

- If the device has not been used for a long time, wipe it down and turn it on to check if it works properly. Then charge it. If any irregularities occur, contact the service.
- Persons with implanted electronic devices, such as pacemakers, must not use the device.
- Using the device near a short wave or microwave processing device (eg 1M) may cause the output of the device to become unstable.
- Do not use this device during pregnancy and menstruation.
- Keep the device and its accessories out of the reach of children.
- Do not immerse the device in any liquid and do not wipe the product with a wet cloth.
- Please turn off the power before charging the device.
- Connect to USB for charging. The device cannot be used while it is charging.
- Improper use or adjustment of the device may void the service warranty.
- Persons with wounds, fractures, scars, cancer, varicose veins and local poor skin elasticity, must not use the device.
- Persons with serious heart, liver, kidney diseases or high fever, must not use the device.

More information about the company and products [page 3 - f7]

- [www.beautifly.eu](http://www.beautifly.eu)
- [www.youtube.com/beautifly.eu](http://www.youtube.com/beautifly.eu)

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**[EN] Rules for handling electrical and electronic equipment:** Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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**[EN] Summary declaration of conformity:** The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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