BEAUTIFLY

Derma Glow

PRO

series

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

f2

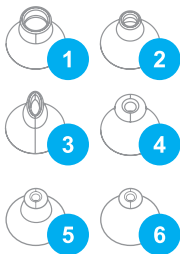f3

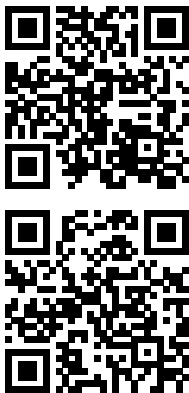Rated voltage:	3.7V
Rated power:	< 1.6W
Battery capacity:	400mAh
Size:	175mm x 46mm x 55.5mm
Weight:	305g
Working time:	> 100 minutes
Charging time:	2 hours
Set includes:	4 interchangeable heads, USB cable, case, replaceable filters, user manual

Safety

Read this manual before using the product and follow the instructions in it when using the device.

Product description [page 2 - f1]

1. Head
2. Display screen
3. Power button
4. USB charging port
5. Back charging indicator

Replaceable suction heads [page 2 - f2]

1. Small round suction head - Facilitates the removal of blackheads from hard to reach places.
2. Round suction head - Vacuum lifting produces a gentle, stimulating massaging effect that improves microcirculation and tightens the skin.
3. Elliptical suction head - It can be used at the eye corner, nose-wing and skins around the mouth.
4. Microcrystalline round suction head - Removes dead and calloused epidermis, acne scars.
5. Microcrystalline tower suction head
6. Microcrystalline small round suction head

Charging

1. Charge the device for 8 hours at most.
2. Make sure the device is off when charging.
3. You can connect the device to a power adapter via USB cable, or directly to a computer.
4. When charging, the back charge indicator flashes, and when fully charged, the charge indicator lights up.

Precautions

- Dry your hands before using the device.
- Do not rinse or immerse the device under water.
- Do the daily skin care after using this device.
- Do not use other devices for mechanical facial cleansing while using this device.
- Do not put too much pressure on your skin while using, as the device itself has a strong suction.
- The suction heads should be disinfected to avoid the emergence and spread of viral, bacterial and fungal diseases of the skin.

Contraindications

- Active purulent acne
- Rosacea
- Retinoid therapy
- Acid treatments, lasers - at least 4 weeks before microdermabrasion
- Psoriasis
- Viral, bacterial and fungal diseases of the skin
- Vascular skin
- Tendency to scarring
- Broken skin (fresh wounds)
- Skin tumors or suspicious moles in the treated area
- Fillers in the face - at least 4 weeks before performing microdermabrasion, otherwise the preparation may move.
- Do not wash the device with water, the device is not waterproof.

More information about the company and products [page 2 - f3]

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EN: Electrical and electronic equipment, as well as batteries and accumulators, must not be mixed with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Used electrical and electronic equipment, as well as batteries and accumulators, should be taken to local waste collection points or returned to the seller, and further details can be obtained from your municipality. Electrical and electronic equipment, as well as batteries and accumulators, may contain hazardous substances, mixtures and harmful components that are dangerous for the environment and human health. Improper handling or damage to these items can lead to health problems or environmental pollution during the disposal or recycling of equipment.

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EN: The manufacturer hereby declares, on their sole responsibility, that the product Derma Glow Pro complies with the Directive of the European Parliament and of the Council 2014/35/EU of 26 February 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, the Directive of the European Parliament and of the Council 2014/30/EU of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and the Directive of the European Parliament and of the Council 2011/65/EU of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU with regard to the list of restricted substances. The full text of the Declaration of Conformity is available at:

<https://beautifly.eu/deklaracje>

Importer: Beautifly Sp. z o.o. Plac Bankowy 2, 00-095 Warszawa, Polska
Model: Derma Glow Pro