BEAUTIFLY

Plasma

Blush

series

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

f1

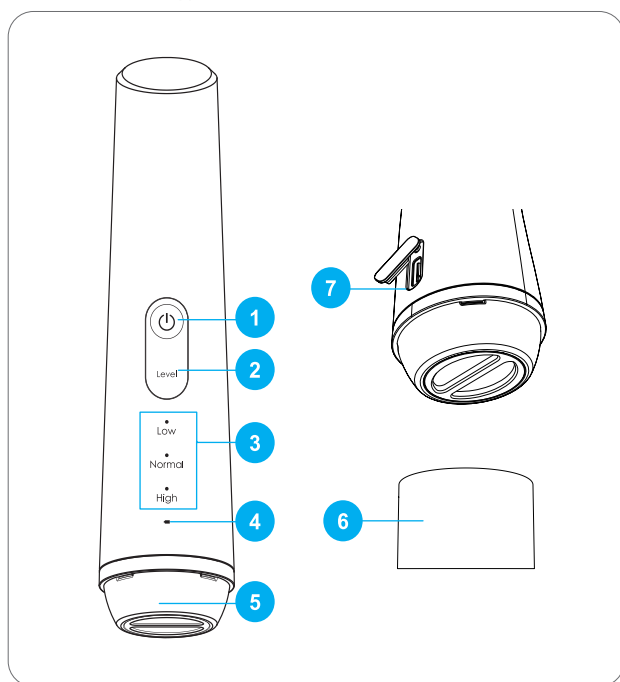f2

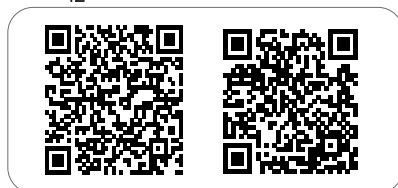Device Size:	163×45×39mm
Net weight:	120 g
Charging way:	charging Wried Type-c
Charging time:	about 3 hours
Battery:	1000 mAh lithium battery
Material:	ABS,PC
Waterproof Standard:	IPX5
Power supply:	AC 100-240V 50/60Hz DC 5V-1A

Security:

Before using the product, read this user manual and follow the recommendations contained therein when using the device.

Product description [page 2-4]:

1. Power on/off button
2. Run button.
3. Gear indicator.
4. Charging indicator.
5. Plasma treatment head.
6. Head cap.
7. Charging port.

Waterproof

The device has an IPX5 waterproof standard, but is not completely waterproof.

Accessories are not waterproof.

Tips for contacting the device with water:

- do not use the product in bathrooms, showers or other humid and wet places
- do not immerse the charging port
- do not immerse or wash the device directly in water
- do not immerse the device in any container filled with water. Keep the device away from damp and wet places

Detailed user manual:

Before using the device, you must charge it. When the charging indicator turns from red to white, charging is complete. After using the product, wipe any remaining moisture from the device.

Charging indicator:

Charging - White light on

Charging complete - Red light flashes

Low power - Charging time

Charging time - about 3 hours

Working time - approximately 60 minutes

Once charging is complete, disconnect the device from the power supply. Leaving the device plugged in after fully charging may reduce battery life.

Power On/Off:

- briefly press the button to turn on the device.
- hold the button for 1 second to turn off the device.
- the device will turn off automatically after 10 minutes.

Changing the intensity level

1. Turn on the device, when turned on, the default intensity level is set to low.
2. Press the <Level> button to increase the intensity level.

Before use:

Before starting the treatment, you should remove make-up from your face and thoroughly clean it of any make-up remnants. Then moisten your face with a cotton pad. During the treatment, the skin must remain moist, so this action must be repeated as often as necessary.

Principles of operation of treatments:

Anti-acne care

Acne-prone skin, in addition to thorough cleansing, requires effective protection against bacteria. Anti-acne therapy helps remove skin irregularities caused by acne: pimples, dimples, subcutaneous nodules and scars. Thus, it restores the uniform structure of the skin and significantly smoothes the complexion. It regulates the keratinization of the sebaceous glands, thereby inhibiting the growth of bacteria and preventing the formation of new skin eruptions. Anti-acne care is the best way to strengthen the protective barrier of the epidermis and eliminate skin inflammation.

Skin regeneration/ intensive rejuvenation

Beautiful, healthy skin at any age? Yes! Regenerating massage using plasma mode slows down the aging process and helps reverse the signs of passing time. Deep revitalization stimulates the cells of each of the three layers of the skin for intensive renewal. It stimulates blood and lymph microcirculation, which rejuvenates even severely damaged skin. Just a few regenerating treatments noticeably reduce wrinkles, providing a delicate lifting effect. Deep regeneration means elastic, firm and velvety smooth skin. Just a youthful look.

Skin care - Plasma mode

Preparation - Before starting the treatment, remove makeup, cleanse the skin thoroughly, and then dry it.

Anti-acne care - Holding the device in your hand, perform circular movements on the area affected by acne. After 10 seconds, perform the treatment on the next area of skin affected by acne. Repeat the treatment 2-3 times for areas requiring intensive care. After completing the treatment, apply a serum adapted to the severity of your acne and gently massage the skin until the preparation is absorbed. No other skin care products are needed.

Skin regeneration/intensive rejuvenation - Hold the device and move it slowly from bottom to top in the direction indicated by the arrow. The movement from the start point to the end point should take approximately 3-5 seconds. The treatment should be repeated 3-5 times. After completing the treatment, apply daily skin care, according to your skin needs and habits. When caring for your skin, do not perform treatments too close to the eye or throat area.

Tips:

1. When using the device for the first time, it is recommended to select a low setting and limit use to 2-3 minutes for the entire face. During subsequent treatments, it is recommended to gradually switch the device to higher settings and extend the duration of use. Do not use the device for longer than 10 minutes on the entire face in one session.
2. The device is characterized by deep absorption (absorption). After treatment, you should choose gentler skin care products (stronger effectiveness of skin care products may cause skin irritation). Be careful with more irritating ingredients such as alcohols, Vc, fruit acids, niacinamide and 377. If you notice any sensitivity or discomfort, discontinue use immediately.
3. It is recommended to use the device at night. If the area affected by acne requires strengthening with treatment during the day, apply sunscreen to the treatment area 1 hour after treatment.
4. There will be sound and smell when using the machine. This is a natural phenomenon that occurs during plasma generation. Continue to use your device freely.
5. During the treatment, the charging port must be closed with a silicone plug. If the cap is securely attached, you may feel a slight current when touched with your hand.
6. If you use the plasma mode, apply the tonic after the treatment.

Acne treatment:

- area of application - acne area
- use time - 5-10 minutes
- frequency - 1-2 times a day

Rejuvenation/penetration into the skin:

- area of use - face
- use time - 5-10 minutes
- frequency - 1 time a day

Contraindications:

- The device should definitely not be used around the eyes, ears, Adam's apple and mouth.
- The device is not intended for children under 15 years of age. Young people aged 15-18 may use the device with the consent and/or supervision of their parents or guardians.
- Do not use the device if the skin in the area to be treated is tanned or red. Using the device on tanned or red skin may cause temporary or even permanent discoloration of the skin.
- Do not use the device if you are pregnant, trying to become pregnant, think you may be pregnant or are breastfeeding. The device has not been tested for this purpose and should therefore not be used in these cases. Hormonal changes that occur during pregnancy or breastfeeding may increase the risk of skin damage or skin sensitivity.
- Do not use the device if you suffer from diseases such as epilepsy, epileptic seizures or migraine. Energy flashes can cause severe headaches and migraines.
- Do not use on areas of skin with tattoos, dark brown or black spots (such as dark freckles, moles or warts). Tattoos or dark spots can absorb too much light energy, which can cause skin swelling, temporary or permanent blisters, discoloration, or skin whitening.
- Do not use the device if you have a tendency to develop keloid scars.
- Do not use if there is any active skin inflammation in the treatment area, including acne, psoriasis, eczema, burns, herpes, abrasions or rash. These factors can make the skin more sensitive to light and make the treatment uncomfortable or cause swelling, blisters or infections on the skin. Once your skin has healed sufficiently, consult your doctor before using the device.
- You should not use the device if you are taking medications that make your skin sensitive to light, including non-steroidal anti-inflammatory agents (e.g. aspirin, ibuprofen, cetaminophen), tetracyclines, phenothiazines, thiazide, diuretics, sulfonamides, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, alpha-hydroxy acids (AHAS), Beta-Hydroxy Acids (BHAs), Retin-A, Accutane and/or topical retinoids.

Precautions

- Do not use the device for purposes other than those for which the product is intended.
- Do not share this product with family members or others. Cross infection may occur.
- Do not use the product in flammable or explosive environments.
- Do not use the device simultaneously with other electrical products (beauty devices, etc.).
- Do not use a power source outside the 100-240V range to charge the device.
- Do not use the product or accessories if they are damaged.
- Do not soak or wash the device directly under water, or immerse it in a sink or bathtub.
- Do not use the product with wet hands.
- Do not drop the main unit of the device, hit it or subject it to strong shocks.
- Do not use a power supply other than the one supplied with the product, as this may cause the product to malfunction.
- If you experience pain, allergies or other adverse reactions during use, stop using the product immediately and consult a doctor.
- Do not place the main unit of the device in places with high temperature, humidity or direct sunlight.

Cleansing:

The device is not waterproof. After each treatment, wipe the product with a dry towel. Do not use corrosive agents to clean the device.

Storage:

The device should be stored in a cool, dry place and out of the reach of children and animals.

More information about the company and products [page 2 - f2]:

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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