BEAUTIFLY

LumiThera epilseries



IPL hair removal
User manual

[f /beautifully.global](#) [/beautifuly_global](#)

[@beautifuly.global](#) [/@beautifuly.global](#)

f1

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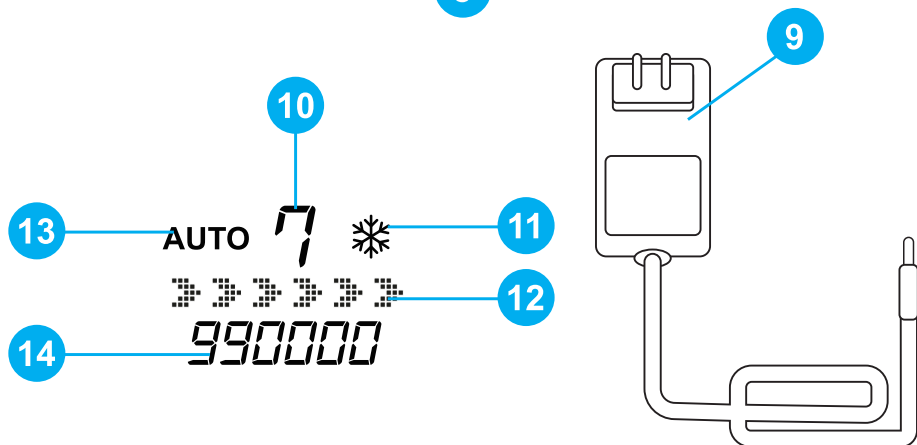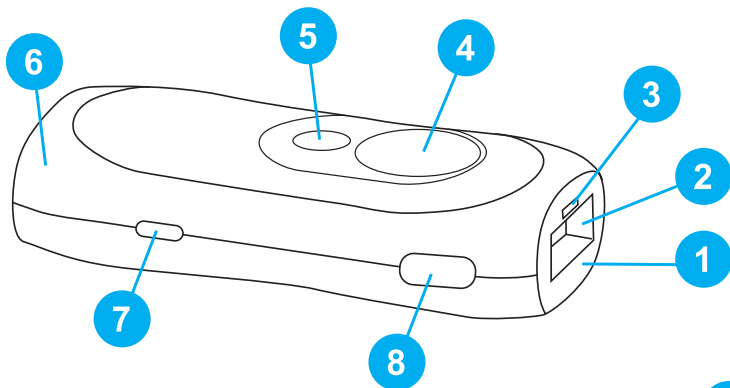Fluence: 1.6 J/cm²~6 J/cm²

Exposure time/pulse duration: 8~12 ms

Interval between flashes only in automatic mode: 1.20~3.02 s

Selectable energy levels: ~9 levels

Exposure area: 5.5 cm²

Included: IPL ×1, user manual ×1, power adapter ×1

Storage and transport conditions

Temperature: -25 degrees ~ 70 degrees

Humidity: 10% to 90% RH

Pressure: 70~106 KPa

Operating conditions

Temperature: 0°C~30°C

Humidity: 10~75% RH

Operating pressure: 70~106 KPa.

Product description (page 2, f1)

1. Cooling surface contacting the skin
2. Flash window
3. Skin color sensor
4. LCD screen
5. Flash button
6. Power socket
7. Power switch/Energy level adjustment switch
8. Air outlet
9. Power adapter
10. Intensity level
11. Cooling symbol
12. Skin detection symbol
13. Automatic flash mode
14. Device lifespan

Contraindications

1. DO NOT use the device if you have very dark skin. Very dark skin may absorb too much light energy and damage (e.g., burn, tattoos, scars, etc.) the skin.
2. DO NOT use the device if your skin has recently been artificially or naturally tanned or is sunburned due to excessive sun exposure.
3. DO NOT use if you are sensitive to sunlight (photosensitivity) or are taking medications causing photosensitivity, e.g., certain antibiotics like tetracycline, NSAIDs like ibuprofen, or retinoids such as Accutane and/or topical retinoids like Retin A.
4. DO NOT use the device if you have a history of skin cancer or precancerous lesions (e.g., moles or a large number of freckles).
5. DO NOT use if you are pregnant, breastfeeding, or menstruating.
6. DO NOT use the device on persons under 18 years of age.

Warnings

1. Do not use the device on the following areas:

- Around the eyes and near the eyebrows,
- On nipples, areolas, labia minora, vagina, anus, and inside nostrils and ears.
- Men cannot use it on the scrotum.
- On and near any artificial elements such as silicone implants, pacemakers, subcutaneous injection ports (e.g., insulin pumps), or dermal fillers.
- On moles, freckles, large veins, darker pigmented areas, scars, and skin lesions without consulting a doctor. This may cause burns and skin discoloration, which can complicate the diagnosis of skin diseases.
- On nipples, tattoos, or permanent makeup; this may cause burns and skin discoloration.
- On areas where you use long-lasting deodorants. This may cause skin reactions.
- Never use this device on sunburned, recently tanned, or artificially tanned

skin. Note: The list is not exhaustive. If unsure whether you can use this device, consult a doctor.

2. Stop using the device immediately if you notice extensive skin redness.
3. Skin damage may occur after prolonged or repeated use of the device in one area.
4. If the power cord is damaged, it should be replaced by the MANUFACTURER, authorized service, or similarly qualified persons to avoid electric shock risk.
5. The device emits intense pulsed light flashes, which can be harmful to the eyes if directly exposed.
6. DO NOT look directly into the flash window when the device is connected to power and turned on.
7. DO NOT try to activate the device towards the eye.
8. It is recommended to wear protective glasses and look away from the device when activating the flash button.
9. The device is intended solely for hair reduction. Do not use it for any other purposes.
10. Using other controls, adjustments, or performing actions not specified in this manual may result in exposure to hazardous radiation.
11. Using an adapter that does not meet specifications is prohibited as it may damage the device.
12. Always check the device before use. Do not use the device or adapter if damaged.
13. Never insert any objects into the device.
14. Do not use the device if the flash window filter is damaged.
15. If the device is damaged, do not touch any internal parts or avoid electric shock.
16. Never use compressed air, rough sponges, abrasives, or aggressive liquids such as gasoline or acetone to clean the device.
17. No part of the device can be serviced or maintained while in use by the patient.
18. Do not disassemble, repair, or modify the product without qualified personnel present, as this may cause damage, malfunction, or injury.
19. Do not leave the device unattended when switched on to avoid fire or burns risk.
20. Unauthorized modifications are prohibited as they may lead to unacceptable risk.
21. Do not use if the device becomes uncomfortably hot to the touch.
22. Do not attempt to activate (flash) the device on any surface other than skin.
23. Do not block the device's air outlet.
24. The device is not intended for use by persons with reduced physical, sensory, or mental capabilities or lack of experience and knowledge unless properly instructed or supervised to avoid fire or burn risk.
25. Do not expose the device to strong shocks to avoid lamp damage.
26. The device is not to be used after use and in case of power failure to avoid device damage.
27. Do not use this device near bathtubs, showers, sinks, or other containers with water.
28. Do not position the device in a way that hinders the operation of the power plug, which is used to disconnect the device from the electrical network.
29. The included power adapter is part of the device; do not use another adapter for power.

Precautions

1. Before connecting the device, check that the voltage on the device matches the mains voltage to avoid electric shock risk or permanent device damage.
2. If the power outlet you use has poor connections, the device plug heats up. Make sure you connect the device to a properly installed outlet to avoid fire and burn risk.
3. If the power adapter is damaged, always replace it with the original type to avoid danger.
4. For hygienic reasons, the device should be used by only one person.
5. When not in use, the device should be stored in a dry place; do not place heavy objects on it.
6. The device and power adapter should be kept in a dry place. The device is not designed for washing.

Indications for use

The IPL hair removal device is an over-the-counter device intended for removing unwanted hair. The device is also indicated for permanent hair reduction, defined as a long-term, stable reduction in the number of regrowing hairs measured at 6, 9, and 12 months after the end of treatment.

Skin tone (page 3, f2)

This device is suitable only for persons with Fitzpatrick skin types I-IV. BEFORE USING THE DEVICE, select your skin tone based on the Fitzpatrick Phototype Scale developed by Harvard Medical School dermatologist Thomas Fitzpatrick, M.D., PhD, used to classify complexion and sunlight tolerance.

The scale has six skin tone levels. Please IDENTIFY your skin tone and check the diagram below to see if the product is suitable.

Function operation instructions

1. Press the power switch for 2 seconds and press again to change the level.
2. Press and hold the flash button for 2 seconds to switch between manual and automatic mode.
3. 990,000 flashes, levels 1-9 consume 1 flash.
4. The display screen first shows 990,000 flashes and decreases with use until zero.

Product functions

IPL hair removal uses intense pulsed light technology which emits multiple wavelengths absorbed by melanin in hair, causing heating. This destroys the hair follicle sufficiently to remove hair and prevent regrowth.

The cooling function of the hair removal device depends on a thermoelectric cooler (TEC). When the device performs hair removal, the TEC creates the Peltier effect. Current flows through a junction formed by two different conductors, electron energy levels jump. The electron must absorb heat, causing the cooling area contacting the skin to create a cold surface, thus lowering the skin temperature in hair removal areas.

Preparation for use

Before use, make sure your skin is in good condition, without injuries or wounds. Skin should be prepared by cleaning, using the included razor to remove all visible hair, and ensuring it is hair-free, completely dry, and not greasy.

Do not use depilatory creams, wax, etc., as chemicals may cause skin reactions. Clean the skin contact surface and filter with a soft cloth. Check if your skin is suitable for this product. (page 3, f2)

Automatic skin tone detection

When powered on, "0" appears on the LCD display; at this moment, place the skin contact surface on the area you want to depilate. If the color symbol on the display stops flashing and "1" or "2" or "3" or "4" or "5" appears, the device is suitable for your skin tone; the number on the LCD is the recommended energy level for your skin tone: level 1 for type V, level 2 for type IV, level 3 for type III, level 4 for type II, level 5 for type I. If the color symbol stops flashing but "0" remains, skin tone detection failed and the device cannot operate. Check the skin tone card above to confirm whether your skin tone is suitable for depilation.

Skin irradiation test

Before first use, WE RECOMMEND testing skin for reactions. Choose a spot on the arm or leg for the test. Connect the power adapter and turn on the device. Set the level to 1, firmly place the skin contact surface on the selected spot and press the flash button to emit light. Move the device to another test spot and increase energy level, press the flash button to perform another test. If the sensation is comfortable, increase the setting by one level and try another test spot. Do not use more than one flash on the same spot. For lighter skin, check results after 30 minutes; for darker skin, after 60 minutes – if no redness occurs, the device is suitable for your

skin. If no skin abnormalities occur but you feel intense heat during the session, the last energy level is most suitable for your skin. (Wear glasses before the test to protect eyes from light.)

Use – steps

1. Before each use, clean the skin contact surface and filter with a soft cloth.
2. Connect the power adapter and press the power button for 2 seconds to turn on the device.
3. When the device is on, "0" appears on the LCD display. You can use the automatic skin tone detection function (see page 3, f2). If the device detects your skin tone is suitable, the recommended energy level appears on the LCD and you can start depilation at this level. You can also manually adjust energy level by pressing the power button to cycle through levels 1-2-3-4-5-6-7-8-9-8-7-6-5-4-3-2-1, as shown on the LCD. It is recommended to start at level one during the first use and gradually increase energy level.
4. Wear professional glasses to avoid eye damage.
5. After setting energy level, place the skin contact surface at a 90° angle to the skin to ensure full contact. If the skin detection symbol on the display stops flashing, press the Flash button to emit light. If the skin contact surface does not fully adhere to the skin, pressing the Flash button will not emit light.
6. After pressing the Flash button, the device emits a single light flash accompanied by a gentle warm tingling. You will feel the heat of the flash. Move the device to the next depilation area. The skin detection symbol on the LCD stops flashing, indicating full skin contact – then repeat the action.

Recommended usage schedule

Weeks 1-4: 3-4 times per week

Weeks 5-6: 2-3 times per week

Weeks 7-8: twice a week
8 weeks and more: 1-2 times per week

Note: The frequency depends on individual hair growth rate. It varies for everyone.

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[EN] Used electrical and electronic equipment, as well as batteries and accumulators, must not be mixed with other waste. This is indicated by the symbol of a crossed-out bin placed on the product packaging. Electrical and electronic equipment, as well as batteries and accumulators, may contain hazardous substances, mixtures, and components harmful to the environment and human health. Improper handling or damage to them can lead to health risks or environmental pollution. Proper handling of used electrical and electronic equipment, as well as batteries and accumulators, in households positively contributes to reuse and recovery, including recycling of the used equipment. Used electrical and electronic devices, as well as batteries and accumulators, can be returned at local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and served the same functions as the equipment sold. When delivering equipment intended for households, the distributor is obliged to collect used equipment from households free of charge at the delivery location, provided that the used equipment is of the same type and served the same functions as the delivered equipment. A distributor operating a retail unit with a sales area, as defined in Article 2, point 19 of the Act of March 27, 2003, on spatial planning and development, of at least 400 m² dedicated to the sale of equipment intended for households, is obliged to accept used equipment from households free of charge at this unit or in its immediate vicinity, provided that none of the external dimensions exceeds 25 cm, without the necessity of purchasing new equipment intended for households. The product was placed on the market after August 13, 2005.

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/UE of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/UE of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/UE of the European Parliament and of the Council as regards the list of restricted substances. The full text of the declaration of conformity is available at the following internet address:
www.beautify.eu/deklaracje

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