BEAUTIFLY

SkinPure Therapy

P R O

series

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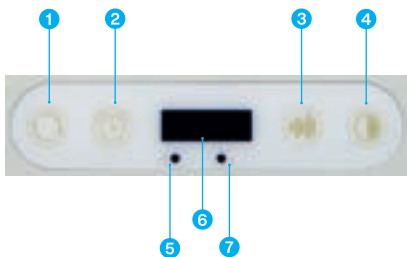
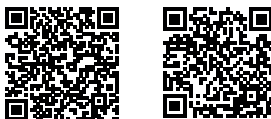
Security:

Before using the product, read this user manual and follow the recommendations contained therein when using the device.

Control panel description [page 2-f1]:

1. Lighting mode selection button: press 1,2,3 times to select red light + infrared light, red light and infrared light. "R I", "R" "I" appears on the display. Then press 4 times to close all lights, the display will show OFF. and so on.
2. Timer function button: press 1,2,3,4 times for 5,10,15,20,25,30 minutes. When "ON" appears on the display, the timer function will not start.
3. Pulse function button: Press 1,2,3 times for 10Hz, 20Hz, 40Hz. The display will show "10", "20", "40". Then press 4 times to close the pulse function, the display will show "on".
4. Dimmer function button: default brightness is 100%: Press 1,2,3,4 times to get brightness 25%, 50%, 75% 100%, the display will show "25", "50", "75", "100".
5. Red light indicator during operation.
6. Digital display showing measurement time and remaining time.
7. Infrared light indicator during operation.

Light sensitivity test:

Before using the device for the first time, check whether you are allergic to light: Place the back of your hand against the device for 3 minutes. If skin becomes red or itchy, stop use and consult a dermatologist. Always clean your face before use to avoid reactions to cosmetics. If baking, stop using the appliance immediately.

LED light therapy

Red and near-infrared light therapy, also known as photobiomodulation (PBM), involves the use of specific wavelengths of natural light in the mid-600 and mid-800 nanometers to improve the health and function of our cells without harmful UV rays, excess heat, or negative effects. Light therapy is one of the earliest registered treatments. Solar therapy was first used by the Egyptians, and forms of light therapy were also practiced by the ancient Greeks, Chinese and Indians. There is no doubt that light has biological effects: in fact, the body needs light to be healthy.

Advantages of light therapy

1. Stimulate hair growth:
LED light therapy stimulates hair growth phases, preventing baldness in men and women.
2. Improve skin health:
LED therapy increases collagen production and tissue repair, improving wound healing, skin appearance and reducing wrinkles.
3. Improve circulation:
NIR light improves circulation, reducing inflammation and relieving pain symptoms in joints and muscles.
4. Reduce recovery time:
NIR accelerates recovery by reducing inflammation and pain, which promotes faster recovery from injuries and illnesses.
5. Reduce Inflammation:
Near infrared (NIR) light lowers inflammatory markers, relieving pain and symptoms of autoimmune diseases and injuries.
6. Promote Healthy Cellular Function:
LED light therapy stimulates collagen production, strengthening hair, skin and connective tissue.

7. Increase Fertility:

LED light can increase testosterone production and improve fertility, positively affecting reproductive function.

8. Pain Reduction:

NIR light therapy effectively reduces neuropathic pain and general discomfort, improving quality of life.

Maintenance and care of the device:

Store this product in a dry place at room temperature and avoid exposure to dust, direct sunlight or humid conditions. Clean the mask before and after use.

Precautions:

Please keep the product away from small children and animals to avoid treating it as a toy. Handle the product carefully when using. This device is designed as a home electronic product and is not suitable for diagnosis or treatment. It should not be used by pregnant women, preparing for pregnancy or breastfeeding, people sensitive to light or taking photosensitizing drugs, with genetically determined skin diseases or conditions such as lupus, photosensitive eczema, or albinism. Do not use while charging. Do not use the device if it is damaged. Do not repair, disassemble or modify the device yourself. Any irregularities should be assessed by the service. Do not drop or hit the device to avoid malfunction. After each treatment, wipe the product with a soft damp cloth. Do not wash the device under water or submerge it. Do not use corrosive agents to clean the device. Store it in a cool, dry place.

More information about the company and products [page 2 - f2]

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[EN] Rules for handling electrical and electronic equipment: Waste electrical and electronic equipment as well as batteries and accumulators cannot be mixed with other waste. This is indicated by the crossed out bin symbol on the product packaging. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage may result in damage to health or environmental pollution. Proper household handling of waste electrical and electronic equipment and batteries and accumulators positively contributes to the reuse and recovery, including recycling, of waste equipment. Waste electrical and electronic equipment as well as batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and performed the same functions as the equipment sold. When delivering equipment intended for households to the buyer, the distributor is obliged to collect used equipment from households free of charge at the place of delivery of the equipment, provided that the used equipment is of the same type and performed the same functions as the delivered equipment. A distributor running a retail unit with a sales area within the meaning of Art. 2 point 19 of the Act of March 27, 2003 on spatial planning and development of at least 400 m² devoted to the sale of equipment intended for households, is obliged to accept waste equipment from households free of charge in this unit or in its immediate vicinity, none of the external dimensions exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market on August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer declares, under its sole responsibility, that the SkinPure Therapy product is in compliance with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, Directive 2014/35/EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment intended for use within specified voltage limits and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances and Directive (EU) 2017/2102 of the European Parliament and of the Council of 15 November 2017 amending Directive 2011/65/ EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment. The full content of the declaration of conformity is available at www.beautifly.eu/deklaracje

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