medivon.

User manual

Collar massager

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1. Configuration and operation


It is recommended to fully charge the device's battery using the charger before use. For this purpose, use only the charger included in the set.

It is recommended not to fully discharge the device to extend the battery life.


It is forbidden to use the massager while charging the battery. This may cause damage or fire.

Controller support:

1. Power button - press to turn on the massager, press again to turn off the massager.
2. Mode Switch Button - Click the button to cycle between 2 modes: Vitality and Relaxation.
3. Heating button - press while the massager is running to turn the heating function on or off.
4. Intensity button - short click the button to switch between 2 massage intensity levels.

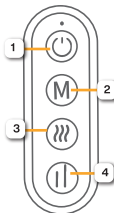2. Using the massager

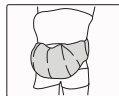
Belly


Neck and shoulder


Thigh

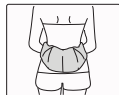Waist

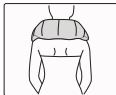Neck

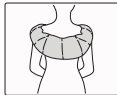Back

3. Safety instructions

Before using the massager, please read the instructions carefully to use the device safely and efficiently. Before first use, please charge the product for approximately 2.5 hours to fully charge it. While the product is charging, it cannot be turned on or used. Before using the massager, please remove jewelry, metal accessories, etc., to avoid discomfort during use. Before using, wait one hour after your last meal.

Keep the massager away from children - it is not a toy! Use the massager for no longer than 10 minutes per session. After this time, the device will turn off automatically. Take a break for a few minutes before you resume work. Unplug the adapter from the socket when not in use.

Do not use the device if you are pregnant or suffer from heart problems or have any physical injuries. Do not use the device if you feel unwell. If you are unsure, consult your doctor before use.

Do not use or leave the massager in the bathroom or other places exposed to moisture. If the device does not work properly, immediately disconnect the power supply and contact the service center. If the power cord is damaged, do not attempt to start the device. The cable must be replaced by the manufacturer, its service agent or a similarly qualified person in order to avoid danger to life or health. Avoid contact with sharp edges or pointed objects that may cut or puncture the surface of the device.

4. Care and cleaning

Always unplug the device before cleaning. Never touch the plug with wet hands. Wait until the device cools down before cleaning.

Do not use wipes containing chemicals or alcoholic or solvent-based liquids. Never immerse any part of the device in liquid.

If your device is particularly dirty, follow these instructions:

1. Dampen a cloth in water or a 3-5% mild detergent solution.
2. Wipe dirty areas with a wet cloth.
3. Rinse the cloth in water and wring it out. Then wipe any detergent residue from the surface of the device (until it dries).

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment. The device was placed on the market after 13.08.2005.

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The manufacturer hereby declares that the Medivon Sooth device complies with Directive 2014/35/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65/EU of the European Parliament and of the Council of June 8, 2011 on restrictions on the use of certain hazardous substances in electrical and electronic equipment. The full content of the declaration of conformity is available at the following internet address: www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska
Model: Medivon Sooth