medivon®

Leg Air Massager

Professional therapy

PL | EN | DE | ES | IT | FR | CZ
SL | PT | LT | RU

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1. Device operation

The device has an air massage function and a heating function. The first function is for the calves and thighs, and the heating function is for the knees. There are three massage modes and three intensity levels of the compression massage function. The massager has a built-in battery - when charging, the intensity indicator will turn blue, green or red. When fully charged, the indicator light will be green. Press the power switch on the top of the controller to enter the standby state.

In standby mode, short press the power button to turn on the device - the indicator of this button will turn red. Repeatedly press the button to change modes: the red LED indicates the mode with heating and the blue LED indicates no heating mode. Hold longer to turn off the device.


While the massager is operating, press the mode button to cycle through the three matching compression massage styles. The fourth mode is the mode without compression massage.

Using the intensity button, we can adjust it to suit our tastes and preferences.

2. Cooling compress

With the device turned off, fully unfold the leg massager. Insert the three pre-chilled (or frozen) gel pads into the pockets located on the inside of the device as shown in the diagram below. Carefully put the device back on your leg to enjoy the cold function. There is no need to turn on the device, however, it is allowed to increase the feeling of cold. It is recommended to test the temperature of the compresses gradually to avoid frostbite at low temperature and skin contact.

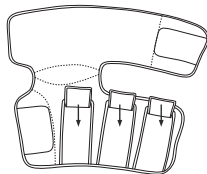3. Dangers

1. Do not use spare parts that are not approved by the manufacturer.
2. Unprofessional service should not disassemble and replace parts of the leg massager without permission.
3. If the power cord is damaged, it must be replaced by appropriate professional personnel to avoid a hazard.
4. Do not use or store the product in a place with high humidity or high temperature.
5. After completing the massage, disconnect the power cord.
6. Do not use the product for more than 15 minutes at a time.
7. If any side effects occur, please stop using the massager immediately.
8. If the temperature is too high or the product starts to smell bad or smells burnt, stop using it immediately and disconnect the power.
9. Do not carry the product by pulling on the power cord to prevent the power cord from loosening.
10. This product is not intended for people who have no physical, sensory or mental abilities or have no experience and knowledge in this field (including children).
11. People insensitive to heat should use the heating function carefully - preferably under the supervision of a caregiver.

4. Warnings

If you have any of the following conditions, consult your healthcare practitioner before use:

- Pregnant women, children or the elderly
- People who suffer from skin infections or have reddened and burned skin
- who use a pacemaker or similar implanted medical device;
- Patients with heart disease or malignant tumors, etc.
- People who are treated or who experience discomfort
- People suffering from osteoporosis
- People with high blood pressure or abnormal blood pressure

5. Cleaning and maintenance

Turn off the product and unplug the power cord before cleaning. Wipe the surface of the product with a dry, soft cloth. Do not clean the product and accessories with harsh detergents such as benzene, thinners, and store in a cool, dry place when not in use. Also avoid exposure to direct sunlight or humid environment, otherwise it may cause color fading or product damage.

Ice pack warning: Store in a well-ventilated place, in a closed container, out of the reach of children. Health hazards: contact with eyes or short hair may cause irritation; ingestion may cause abdominal pain, discomfort. Precautions: Wash thoroughly after handling, do not eat, drink or smoke when handling this product. Wear protective gloves/protective clothing/eye protection/face protection. In case of: contact with skin or hair, immediately remove all contaminated clothing and wash it; rinse the skin with water; in eyes: rinse carefully with water for several minutes; ingestion: rinse mouth, call for medical help.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment.

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The manufacturer declares under his sole responsibility that the Noble Flex device complies with Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, as well as Directive 2011 of the European Parliament and of the Council /65/EU of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to the Directive of the European Parliament and 2011/65/EU as regards the list of restricted substances. The full text of the declaration of conformity is available at the internet address

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska
Model: Noble Flex