BEAUTIFLY

Looky
Young
series


Eye massager User manual

 /beautifly.eu  @beautifly.eu  beautifly.eu

 /beautiflyeu  @beautifly.eu


f1

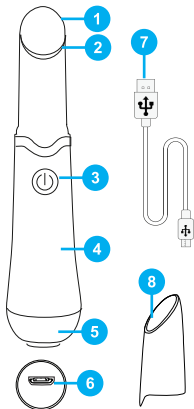f2

f3

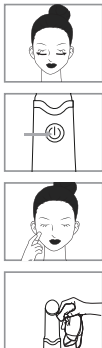f5

f4

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- www.beautify.eu

- www.youtube.com/beautify.eu

Dimensions:	25*25*120mm
Net weight:	30g
Temperature:	42 ± 2°C
Power:	≤3W
Charging voltage:	≡5V
Voltage:	≡3.7V
Charging current:	160mA
Battery capacity:	200mA

Safety: Before using the product, please read this manual and follow the instructions contained in it during use.

Device description [page 2 - f1]

1. Massage head / 2. Red-blue LED ring / 3. Power button/mode switch / 4. Device body / 5. Induction field / 6. USB charging interface / 7. USB charging cable / 8. Protective cover

User manual [page 2 - f3]

1. Apply an appropriate amount of care product to the skin around the eyes. 2. Press and hold the power button for 2 seconds to turn on/off the device. Press the power button once to select the red or blue light therapy mode, depending on your skin care needs.
3. Hold the device so that the massage head is in constant contact with the skin and vibrates gently.
4. After use, clean the massage head with a dry tissue and put on the protective cover.

Note: The device automatically turns off after 2 minutes of non-use

How to use [page 2 - f4]

1. Eye care: Apply eye cream and slowly massage the skin from the inner to the outer corner of the eye. This will improve blood and lymph circulation and remove dark circles. The massage will help the eye cream to be absorbed deeper, providing better hydration and effect on the skin.
2. Forehead care: Apply skin care product and massage in an upward motion - from the eyebrow to the hairline. Ion therapy and high-frequency vibration deeply activate skin cell tissues, improving skin tension above the nose, between the eyebrows and on the forehead.
3. Cheek care: Apply care product and massage the skin from the inside of the face outward to improve blood circulation, achieve a skin rejuvenation effect and minimize the appearance of fine lines.
4. Lip care: Apply lip balm along the lip contour. Gently massage the corners of the lips upwards to eliminate fine lines around them and improve moisture, making the lips soft and shiny.
5. Neck care: Apply moisturizer and massage from the bottom to the top. Ion therapy and high-frequency vibration massage care for the skin, used to revitalize cells, deeply nourish the skin and soften neck lines.

Storage

1. Store the device in a dry place, protected from moisture and erosive fluids.
2. Keep the device out of reach of children.
3. Avoid massaging the eyeballs and do not use the device on injured, bleeding, irritated or swollen skin.
4. Protect the device from falling, strong impacts and sunlight.
5. The device is not waterproof, do not immerse it in water.
6. Keep the device away from damp places and strong magnetic fields.

7. Do not attempt to repair the device yourself - repairs should only be performed by qualified specialists.

8. Do not use the device while charging. The red light flashes when the battery is low and while charging, and a solid blue light indicates a full charge.

Precautions: Before using the purchased product, read the product user manual and these safety warnings, maintenance and storage instructions, and the rules for handling used electrical and electronic equipment. Follow the recommendations contained therein and use the product in accordance with its intended use. The product is not a toy and is not intended for use by children. The product may only be used by adults. Do not use the product if it does not work or works incorrectly, or has been dropped or damaged in any other way. In such a situation, do not open or repair the product yourself. If the product is still subject to liability for defects in the sold item (warranty) or warranty, contact the seller or guarantor. If the warranty or guarantee period has expired, or the defect is not caused by a cause inherent in the product, return the product to the appropriate service for repair. Do not use the product near sources of heat, water or moisture, and do not immerse it in water or other liquids. The product is intended for use at room temperature, in the range of 10 to 40 degrees Celsius. Do not expose the product to direct sunlight. After use, store the product in a dry place, at room temperature, between 10 and 40 degrees, out of reach of children and animals. If the product becomes dirty, disconnect it from the power supply, wipe it with a damp cloth and leave it to dry before storing it. Do not use abrasive materials or corrosive cleaning agents, as this may damage the product. Do not use the product if you notice a damaged power cord or plug, if it emits an unpleasant smell (e.g. burning), overheats or smoke. In such a situation, do not open or repair the product yourself. If the product is still covered by liability for defects in the sold item (warranty) or guarantee, contact the seller or guarantor. If the warranty or guarantee period has expired, or the defect is not caused by a cause inherent in the product, return the product to the appropriate service for repair. Do not leave the product unattended while it is switched on. The power source must match the required voltage of the product. After use, disconnect the product from the power source and allow it to cool down before moving or storing it. Do not twist or coil the power cord, as this may damage it. Accessories for the product that do not come from the manufacturer may cause the product to malfunction or be damaged.

More information about the company and products [page 2 - f5]

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances. The full declaration of conformity is available at <https://beautifly.eu/deklaracje>.

Importer: Beautifly Sp. z o.o., Plac Bankowy 2, 00-095 Warszawa, Polska
Model: Looky Young
Made in China.

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[EN] Rules for handling electrical and electronic equipment: Used electrical and electronic equipment, as well as batteries and accumulators, must not be disposed of together with other waste. This is indicated by the crossed-out bin symbol placed on the product packaging. Electrical and electronic equipment, as well as batteries and accumulators, may contain hazardous substances, mixtures, and components that are harmful to the environment and human health. Improper handling or damage may result in health risks or environmental pollution.

Proper handling of used electrical and electronic equipment, as well as batteries and accumulators, in households contributes positively to reuse and recovery, including recycling, of used equipment.

Used electrical and electronic equipment, as well as batteries and accumulators, can be returned at local selective waste collection points. The distributor is obliged to accept used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and served the same functions as the equipment being sold. When delivering household equipment to the buyer, the distributor is obliged to collect used equipment from households free of charge at the delivery location, provided that the used equipment is of the same type and served the same functions as the delivered equipment. A distributor operating a retail unit with a sales area, as defined in Article 2(19) of the Act of 27 March 2003 on spatial planning and development, of at least 400 m² dedicated to the sale of household equipment, is obliged to accept used equipment from households free of charge at that unit or in its immediate vicinity, provided that none of the external dimensions of the used equipment exceed 25 cm, without the requirement to purchase new household equipment.

The product was placed on the market after August 13, 2005.