BEAUTIFLY

Plasma
PRO
series

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

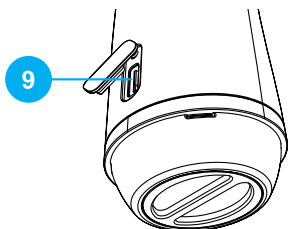
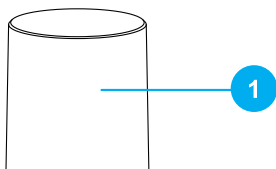
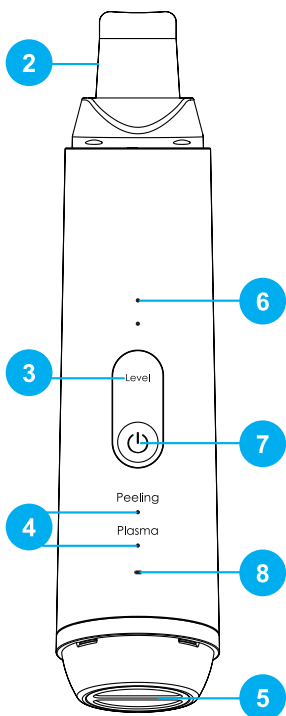
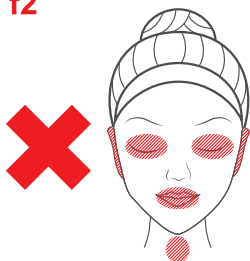
f3

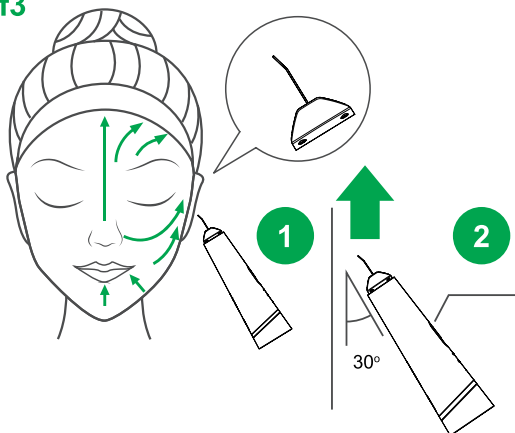f4

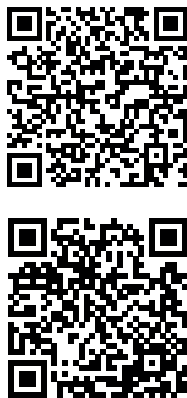Device Size:	175,5x43,5x22mm
Net weight:	141g
Charging way:	USB-C charging
Charging time:	about 3 hours
Battery:	1000 mAh lithium battery
Material:	ABS, stainless steel, PC
Waterproof Standard:	IPX5
Power supply:	AC 100-240V 50/60Hz DC 5V-1A

Security

Before using the device, read the entire Safety Guide. Please follow this User Manual when using the product.

Product description [see page 2, f1]

1. Spatula plug
2. Spatula
3. Run button
4. Operating mode indicator
5. Plasma treatment head
6. Gear indicator
7. Powermode switch button
8. Charging indicator
9. Charging port

Waterproof

The device has an IPX5 waterproof standard, but is not completely waterproof.

Accessories are not waterproof.

Tips for contacting the device with water:

- Do not use the product in bathrooms, showers or other damp or wet places.
- Do not immerse the charging port.
- Do not immerse or wash the device directly in water.
- Do not immerse the device in any container filled with water. Keep the device away from damp and wet places.

Landing

Before using the device, you must charge it. When the charging indicator turns from red to white, charging is complete. After using the product, wipe any remaining moisture from the device.

Charging indicator

Loading > Red light flashes
Charging complete > White light on
Low power > Red light flashes

Charging time > Approximately 3 hours

Working time > Approximately 60 minutes

Once charging is complete, disconnect the device from the power supply.

Leaving the device plugged in after fully charging may reduce battery life.

Power on/off

- Briefly press the <power> button to turn on the device.
- Hold the button for 1 second to turn off the device.
- In standby mode, if no other operations are performed within 15 seconds, the device will turn off automatically.
- The device will turn off automatically after 10 minutes of treatments in plasma and peeling mode.

Function switch

Press the power/function button to select different operating modes: Standby > Peeling/L > Plasma (switching loops)

Changing the intensity level

Press the intensity level button to select the appropriate intensity level.

Peeling mode > 1 indicator light on > Low level

Peeling mode > 2 indicator light on > High level

Plasma mode > 1 indicator light on > Low level

Peeling mode > 2 indicator light on > High level

Before use:

Before starting the treatment, you should remove make-up from your face and thoroughly clean it of any make-up remnants. Then moisten your face with a cotton pad. During the treatment, the skin must remain moist, so this action must be repeated as often as necessary.

Contraindications for use in certain places of the body. [see page 2, f2]

You should definitely not use the device around the eyes, ears, Adam's apple or mouth, as shown in the shaded part of the figure on the right.

Anti-acne care

Acne-prone skin, in addition to thorough cleansing, requires effective protection against bacteria. Anti-acne therapy helps remove skin irregularities caused by acne: pimples, dimples, subcutaneous nodules and scars. Thus, it restores the uniform structure of the skin and significantly smoothes the complexion. It regulates the keratinization of the sebaceous glands, thereby inhibiting the growth of bacteria and preventing the formation of new skin eruptions. Anti-acne care is the best way to strengthen the protective barrier of the epidermis and eliminate skin inflammation.

Skin regeneration/Intensive rejuvenation

Beautiful, healthy skin at any age? Yes! Regenerating massage using plasma mode slows down the aging process and helps reverse the signs of passing time. Deep revitalization stimulates the cells of each of the three layers of the skin for intensive renewal. It stimulates blood and lymph microcirculation, which rejuvenates even severely damaged skin. Just a few regenerating treatments noticeably reduce wrinkles, providing a delicate lifting effect. Deep regeneration means elastic, firm and velvety smooth skin. Just a youthful look.

Peeling

The basic step in proper skin care is to cleanse the skin of dead epidermis cells. Exfoliation by peeling gives excellent and quick results. Visibly improves the appearance of the skin, smoothes and brightens it. Peeling ensures proper keratinization, i.e. keratinization and renewal of skin cells. Thanks to this, it eliminates the risk of developing bacterial skin diseases. Peeling treatments help remove blackheads and additionally facilitate the absorption of nutrients into the skin. The result of the treatments is clean, well-groomed and smooth skin.

Peeling mode [see page 2, f3]

- Before starting the treatment, moisten the skin with water.
- Do not press the spatula hard against the skin or use it repeatedly on the same area.
- If the skin begins to dry out during the treatment, it should be moisturized again with water.

Tilt the back of the skin paddle at an angle of 30° and move it upwards in the direction indicated by the arrow.

Pay attention to the direction of movement of the blade [see page 2, f3]

1. Skin surface
2. Power button outside

Tips:

- It is best to use the device when the skin is warm: after a bath or after putting a hot towel on your face.
- After using the device, take care of your skin by using a tonic.
- If you use the plasma mode, apply the tonic after the treatment.

Anti-acne care

Holding the device in your hand, make circular movements on the area affected by acne. After 10 seconds, perform the treatment on the next area of skin affected by acne. Repeat the treatment 2-3 times for areas requiring intensive care. After completing the treatment, apply a serum adapted to the severity of your acne and gently massage the skin until the preparation is absorbed. No other skin care products are needed.

Skin regeneration/Intensive rejuvenation

While holding the device, move it slowly from bottom to top in the direction of the arrow. The movement from the start point to the end point should take approximately 3-5 seconds. The treatment should be repeated 3-5 times. After completing the treatment, apply daily skin care, according to your skin needs and habits.

When caring for your skin, do not perform treatments too close to the eye or throat area.

Tips:

1. When using the device for the first time, it is recommended to select a low setting and limit use to 2-3 minutes for the entire face. During subsequent treatments, it is recommended to gradually switch the device to higher settings and extend the duration of use. Do not use the device for longer than 10 minutes on the entire face in one session.
2. The device is characterized by deep absorption (absorption). After treatment, you should choose gentler skin care products (stronger effectiveness of skin care products may cause skin irritation). Be careful with more irritating ingredients such as alcohols, Vc, fruit acids, niacinamide and 377. If you notice any sensitivity or discomfort, discontinue use immediately.
3. It is recommended to use the device at night. If the area affected by acne requires strengthening with treatment during the day, apply sunscreen to the treatment area 1 hour after treatment.
4. There will be sound and smell when using the machine. This is a natural phenomenon that occurs during plasma generation. Continue to use your device freely.
5. During the treatment, the charging port must be closed with a silicone plug. If the cap is securely attached, you may feel a slight current when touched with your hand.
6. If you use the plasma mode, apply the tonic after the treatment.

Acne treatment > Application: acne areas > Usage time 5-10 minutes > Frequency: 1-2 times a day

Rejuvenation / Penetration deep into the skin > Application: face > Usage time 5-10 minutes > Frequency: 1 time a day

Contraindications

- The device is not intended for children under 15 years of age. Young people aged 15-18 may use the device with the consent and/or supervision of their parents or guardians.
- Do not use the device if the skin in the area to be treated is tanned or red. Using the device on tanned or red skin may cause temporary or even permanent discoloration of the skin.
- Do not use the device if you are pregnant, trying to become pregnant, think you may be pregnant or are breastfeeding. The device has not been tested for this purpose and should therefore not be used in these cases. Hormonal changes that occur during pregnancy or breastfeeding may increase the risk of skin damage or skin sensitivity.
- Do not use the device if you suffer from diseases such as epilepsy, epileptic seizures or migraine. Energy flashes can cause severe headaches and migraines.
- Do not use on areas of skin with tattoos, dark brown or black spots (such as dark freckles, moles or warts). Tattoos or dark spots can absorb too much light energy, which can cause skin swelling, temporary or permanent blisters, discoloration, or skin whitening.
- Do not use the device if you have a tendency to develop keloid scars.
- Do not use if the treatment area has any active skin inflammation, including acne, psoriasis, eczema, burns, herpes, abrasions or rash. These factors can make the skin more sensitive to light and make the treatment uncomfortable or cause swelling, blisters or infections on the skin. Once your skin has healed sufficiently, consult your doctor before using the device.
- Do not use the device if you are taking medications that make your skin sensitive to light, including non-steroidal anti-inflammatory agents (e.g. aspirin, ibuprofen, cetaminophen), tetracyclines, phenothiazines, thiazide, diuretics, sulfonyleureas, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, alpha-hydroxy acids (AHAS), Beta-Hydroxy Acids (BHAS), Retin-A, Accutane and/or topical retinoids.

Precautions

- Do not use the device for purposes other than those for which the product is intended.
- Do not share this product with family members or others. Cross infection may occur.
- Do not use the product in a flammable or explosive environment.
- Do not use the device simultaneously with other electrical products (beauty devices, etc.).
- Do not use a power source outside the 100-240V range to charge the device.
- Do not use the product or accessories if they are damaged.
- Do not soak or wash the device directly under water, or immerse it in a sink or bathtub.
- Do not use the product with wet hands.
- Do not drop, hit, or subject the unit to strong shocks.
- Do not use any power adapter other than the one supplied with the product, as it may cause the product to malfunction.
- If you experience pain, allergies or other adverse reactions during use, stop using the product immediately and consult your doctor.
- Do not place the main unit of the device in places with high temperature, humidity or direct sunlight.

Cleaning and maintenance

1. Protect the device against moisture and high temperatures, do not place it in the dishwasher or microwave.
2. Turn off the device before cleaning.
3. Do not use harsh, corrosive or volatile cleaning materials.
4. Use a soft cloth or cotton pad moistened with cold or warm water to gently wipe the device.
5. Do not rinse the device under running water.
6. If the device is not in use for a long time (approximately 3 months), remove and fully charge before further storage. This will avoid damaging the battery.
7. Store in a dry, shaded place out of reach of children and animals.

More information about the company and products [page 2 - f2]

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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