BEAUTIFLY

Booster

PRO

series

Fale radiowe
Rejuvenation radio
frequency

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/beautifly.eu

@beautifly.eu

beautifly.eu


/beautiflyeu

@beautifly.eu

f1

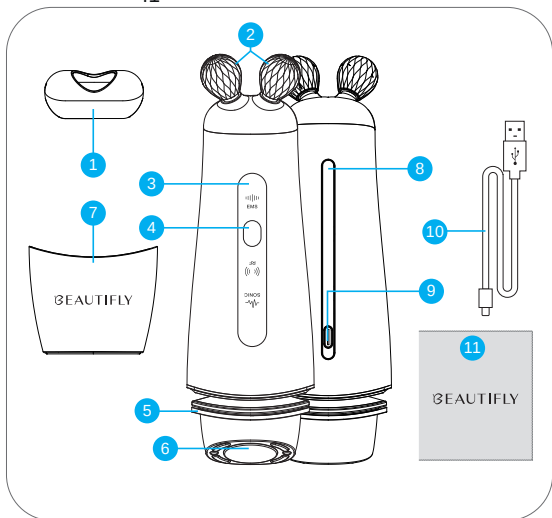f2

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Nominal Voltage:	3.7V
Working Current:	SONIC 100-400mA / RF 200-2800mA / EMS 8-30mA
Ultrasound Frequency:	3MHz±10%
Radio Waves:	1MHz
Water Resistance:	IPX5
Time Modes:	3Min, 5Min, 6Min
Working Time:	<3H
Battery Capacity:	1000 mAh

Description of parts [see page 2, f1]

- | | |
|---------------------------------|-------------------------|
| 1. Transparent protective cover | 7. Magnetic base |
| 2. V-contour roller | 8. Inductive strip |
| 3. HD Display | 9. Type C charging port |
| 4. Mode switching button | 10. Type C cable |
| 5. Magnetic ring | 11. Cleaning cloth |
| 6. Main head | |

SONIC Mode

Press and hold the mode switching button for 1 second, an audio signal will be generated once, the corresponding indicator will light up and automatically turn off after 3 minutes.

RF mode

In SONIC mode, short press the mode switch button to switch, the buzzer will beep once and the corresponding indicator will be on for a long time and will turn off automatically after 5 minutes.

EMS mode

In RF mode, short press the mode switch button to switch, the buzzer will beep once and the corresponding indicator will be on for a long time and turn off automatically after 6 minutes.

Turning off

Manual Power Off: In any mode, press and hold the mode switch button for 1 second to power off. Automatic shutdown will occur after the device has finished operating in any mode and will turn off after an acoustic signal.

Not enough power

When the battery level is less than 10%, the white battery indicator will flash slowly once every 3 seconds and cyclically 3 times. The white slow flash will change to fast flash for the third time and the buzzer will beep twice until it turns off.

In RF import mode, you can use the device 3 times without charging. The low battery indication is the same as the previous two modes.

When the battery level is 0%, the buzzer will beep 3 times, the battery will be turned off and the white indicator will flash for 5 seconds. Then the power will be turned off. The device cannot start.

Battery charging

The white battery light flashes while charging.

Battery charged

The white battery indicator is on.

Memory function

After a manual or scheduled shutdown, the default is the last step mode on the next startup.

How to use

Below is a demonstration of the process (figure and text below). You can choose any mode to use it alone or follow all the steps one by one according to the introduction of "Function Description".

SONIC mode [see page 3, f3]

1. Before use, remove your makeup and cleanse your face with a facial cleanser.
2. Firmly clamp a cotton swab between the magnetic ring and the main head and moisten it with micellar liquid.
3. Turn on the power and enter SONIC mode.
4. Point the main head at your facial skin and move each part of the face and neck from outside to inside to thoroughly cleanse it.
5. On the forehead, nose, lips, cheeks and other parts more exposed to pollution, you can repeat the procedure several times.
6. After 3 minutes, the device will turn off automatically.
7. Remove the cotton pad and clean the main head with a cloth.
8. Rinse your face with clean water and apply tonic to your skin.

RF mode [see page 3, f4]

1. Apply the care essence, serum, ampoule or mask evenly to the face and neck.
2. Power on and switch to RF mode.
3. Slowly move the main head along the contour of the face and apply the treatment evenly to the skin (do not stay in the same position on the skin for too long).
4. After 5 minutes, the device will turn off automatically.
5. Use a cloth to clean the main head.

EMS mode [see page 4, f5]

1. Spread firming care products evenly over your face. It can be a cream that we use to complete the care.
2. Turn on and go to "EMS" mode.
3. Adjust the V face roller to the contour of your face, massage and lift evenly from bottom to top.
4. After 6 minutes, the product will turn off automatically.
5. Use a cleaning cloth to clean the face roller V.

Contraindications

1. Persons under 18 years of age cannot use the Booster device.
2. It is prohibited to use the device for electronic devices inside the body. e.g. pacemaker.
3. People who have undergone any surgery or surgery in the area where the procedure is to be performed in the last month. It is not recommended to use the device after skin exposure to the sun.
4. People with damaged skin, suppuration, inflammation, fresh scars and electromagnetic sensitivity should not use the B-Booster device.
5. People suffering from blood diseases, tuberculosis, celiac disease, malignant tumors, hypertension, fever are not allowed to use the device.
6. Do not use if pregnant, menstruating or breastfeeding.
7. Before starting the treatment, remove watches, glasses, jewelry and other items containing metal.
8. Do not use around the eyeball, mouth or throat.
9. Do not use on the same area of skin for too long.
10. Do not leave this product unattended while it is turned on.
11. Do not use in high temperature environment.
12. If the product is hit or dropped on the ground, check for damage before further use.
13. Check that the plug and power cord are intact to avoid the risk of electric shock.
14. Do not disassemble this device yourself.
15. Please do not turn on the device while charging.

Cleaning methods

After each use, rinse the main head and V face rollers with clean water and then dry with a clean cloth. Do not use alcohol or corrosive liquids for wiping the magnetic ring.

More information about the company and products [page 2 - f2]: www.beautifyli.eu;
www.youtube.com/beautifyli.eu

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under his sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of legislation. Member States in the field of electromagnetic compatibility and Directive 2011/65/EC of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment and Commission Delegated Directive (EU) 2015/863 of 2011. On 31 March 2015, amends Annex II to Directive 2011/65/EU of the European Parliament and of the Council on the list of restricted substances.

Model: Booster Pro