| EN | NB | FI | SV | DA |

BEAUTIFLY

B – Bubble

B o d y
series

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
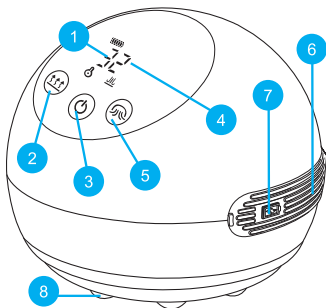
f4

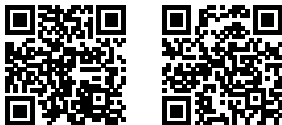Technical Specifications

- Dimensions: 80 x 75 mm
- Weight: 235 g
- Voltage: 3.7V
- Current: 2 A
- Power: 4 W
- Red light temperature: 36-45°C
- Number of gears: 5

Safety rules

Before using the device for the first time, please read the user manual carefully and follow its recommendations to ensure safe use.

Device components (page 2 - f1)

1. Operating mode indicator with temperature information
2. Temperature adjustable
3. On/off button
4. Vacuum mode indicator
5. Vacuum adjustment option
6. Anti-slip silicone surface
7. Charging port
8. Eight massage balls

Red light therapy

The treatment using red light with a wavelength of $625 \text{ nm} \pm 3 \text{ nm}$ effectively combats skin laxity and improves its density. Stimulates the production of collagen and elastin, increasing cell viability and improving blood circulation.

Vacuum massage

Massage using negative pressure stimulates the deep layers of the skin and muscles, supporting regeneration, improving circulation and reducing fat tissue.

Magnetotherapy

The use of an energy magnet supports the body's magnetic balance, purifying the blood and stimulating flow in internal organs.

How to use (page 2 - f2)

Before starting the treatment, apply care oil to the skin. To prevent oil from entering the device, before use, insert a cotton filter into the filter cap and then place it in the suction hole.

User manual (page 2 - f3)

1. To turn on the device, press and hold the button [figure 1]. The device will enter standby mode, which will be indicated by the display 00 on the screen.
2. The device is also turned off by long pressing the button [figure 1].

3. Briefly press the button [figure 2] to activate the vacuum mode (the icon [figure 3] will appear on the screen). Short presses of this button allow you to change the vacuum gears (0-1-2-3-4-5-0 off-1 recycle). The vacuum turns off when the gear is set to 0.
4. To activate red light heating, briefly press the button [figure 4]. The temperature icon [Figure 5] will appear on the screen and the red light will glow. Short presses of the button [figure 4] change the heating mode (0-1-2-3-4-5-0 shutdown-1 recycle). The heating turns off when the device is in gear 0.
5. a. In standby mode, the device will automatically turn off after 5 minutes of inactivity.
 - b. Vacuum mode allows you to adjust the intensity from 1st to 5th gear.
 - c. In temperature mode, it is possible to adjust the heating intensity from 1st to 5th gear. The temperatures of individual gears are: $36^{\circ}\text{C} \pm 2$ (gear 1), $38^{\circ}\text{C} \pm 2$ (gear 2), $40^{\circ}\text{C} \pm 2$ (gear 3), $42^{\circ}\text{C} \pm 2$ (gear 4), $45^{\circ}\text{C} \pm 2$ (gear 5).

Charging your device

After connecting the device to a power source via a USB port, the device stops being active and starts charging. Do not use the device while charging.

The device is equipped with a lithium battery. Please use DC3.7V, 2A power adapter or computer USB port to charge.

NOTE: The device does not work while charging.

Post-treatment tips

After the treatment, you should drink two glasses of water to replenish fluids and support your metabolism. In case of discomfort after the procedure, discontinue use for a few days.

Cleaning and maintenance

Clean the device with a damp towel, avoiding liquids getting inside. Medical alcohol with a concentration of 75% can be used for disinfection.

Storage

Store the device in a dry, cool place, away from children and direct sunlight. Charge the battery regularly if the device is not used for a long time.

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled.

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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