BEAUTIFLY

Scrub Peeling

Young
series

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

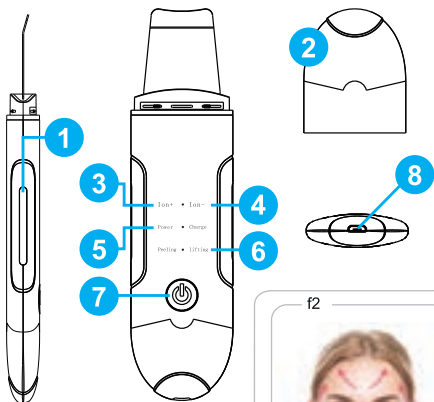
f3

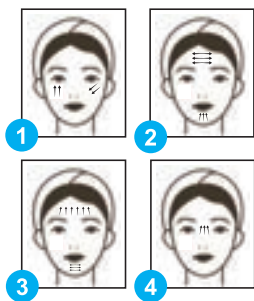f4

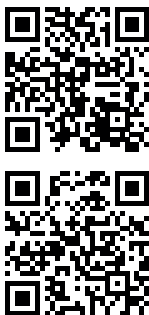Size:	175*43*14mm
Libra:	95g
Ultrasounds:	25KHZ
Material:	ABS
Rated current:	0.65A
Rated voltage:	3.7V
Rated power:	1.5W
Battery Capacity:	500mAh
Included:	device, USB cable, user manual

Security

Please read these instructions before using the product and follow the recommendations contained therein when using the device.

Device description [page 2 - f1]

1. Ion sensor
2. Protective cap
3. Ion+
4. Ion-
5. EMS indicator
6. Ultrasound indicator
7. Power button/mode switch
8. USB charging port

Description of modes [page 2 - f2]

Ion+ ion export mode

Removes impurities from deep pores using positive ions

Ion - negative ion input mode

It introduces nutrients contained in cosmetics deep into the skin.

Power Charge - EMS pulse mode

It lifts and tightens the skin, stimulates dormant collagen tissues, giving a rejuvenating effect.

Lifting peeling - ultrasonic cleansing mode

High-frequency vibrations cleanse the lower layer of the skin, soften and remove dead skin.

Cleansing zones [page 2 - f3]

1. Cheeks
2. / 3. Forehead and jaw
4. Nasal bridge

User manual

1. Hold the power button for 2 seconds to turn on the device. The default setting is lifting mode and exfoliation and the LED1 indicator is always on.
2. Briefly press the power button to cycle through 4 operating modes: lifting mode and exfoliation, LED1 lights up → power charging mode, LED2 flashes quickly → ion import mode+, LED3 lights up → ion import mode-, LED4 lights up.
3. Hold the power button for 2 seconds to turn off the device.
4. The buzzer will beep once when the device is turned on or the operating mode changes. The same sound will appear twice when

the device is turned off.

Landing

1. While charging, LED1, LED2, LED3, LED4 flash in rest mode.
2. When the device is fully charged, LED1, LED2, LED3, LED4 turn off. When the battery level is low, the buzzer beeps twice and the light turns off.

How to use

Epidermis cleansing (Ultrasounds)

1. Before starting the treatment, clean your skin thoroughly and remove makeup.
2. Then moisten your facial skin with a cotton pad or atomizer. For this purpose, you can use physiological saline, tonic, water or hydrolate. During the treatment, the skin must remain moist at all times. The treatment cannot be performed on dry skin.
3. Turn on the device and activate SONIC mode.
4. Gently move the spatula from bottom to top, along the contours of the face at a 45-degree angle. Repeat the action thoroughly several times over the entire face.

Epidermis cleansing with ionization (Ultrasound + Ionization +)

1. Before starting the treatment, remember to keep the skin wet at all times. For this purpose, you can use physiological saline, tonic, water or hydrolate. The treatment cannot be performed on dry skin.
2. Turn on the device and activate ION+ mode.
3. Gently move the spatula from bottom to top, along the contours of the face at a 45-degree angle. Repeat the action thoroughly several times over the entire face.
4. Wipe your face with tonic.

Introducing nutrients using ultrasound (Ultrasound + Ionization-)

1. After cleansing your facial skin, apply a care product (cream, serum, ampoule, gel).
2. Turn on the device and activate ION- mode.
3. Gently move the spatula from bottom to top, along the contours of the face. The spatula should be placed in the opposite direction in step one. Repeat the action thoroughly several times over the entire face.

Note: If you have any doubts about how to hold the device correctly, please visit our YouTube channel [/beautify.eu](https://www.youtube.com/channel/UCbeautify)

EMS (Ultrasound + EMS Currents).

1. After caring for your facial skin, turn on the device and activate EMS mode.
2. To improve skin elasticity, move the spatula from bottom to top, along the contours of the face at a 45-degree angle. Repeat the action thoroughly several times over the entire face.
3. The device will turn off automatically after 10 minutes of operation
4. After finishing using the device, wipe it with a cotton cloth a damp cloth.

Tips

1. Hold the conductive strip with your hand during use on the side of the device so that the ion import and export function can function properly.
2. When in use, the spatula should work at an angle of approximately 45 degrees as shown in the picture.
3. Do not stop the device during use in one place, but move it slowly.
4. If inflammation occurs on the face, wait until the symptoms subside before use.
5. It is recommended to use the device 2-3 times a week. It should not be used several times in one day.
6. The spatula and body should be clean during operation and secured.

Cleaning and maintenance

1. Protect the device from moisture, do not put it in the dishwasher or microwave.
2. Turn off the device before cleaning.
3. Do not use harsh, corrosive or volatile cleaning materials.
4. Use a soft cloth or cotton pad moistened with cold or warm water to gently wipe the device.
5. Do not rinse the device under running water.
6. If the device is not in use for a long time (approximately 3 months), remove and fully charge before further storage. This will avoid damaging the battery.

More information about the company and products [page 2 - f4]

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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Importer: Beautifly Sp. z o.o.

Plac Bankowy 2, 00-095 Warszawa, Polska

Model: Scrub Peeling Young