# BEAUTIFLY

Booster  
B o d y  
series

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

f1

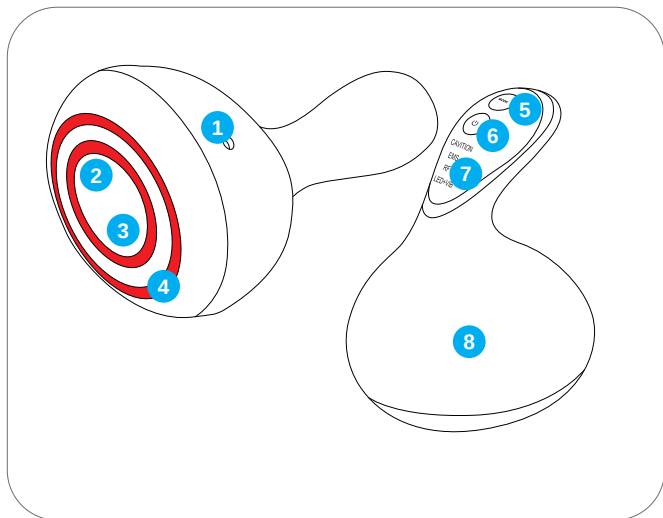f2
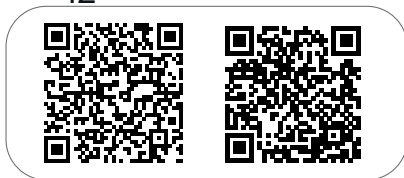
|                  |          |
|------------------|----------|
| Entry:           | 110-240V |
| Exit:            | 5V 1A    |
| Lithium battery: | 900mAh   |
| Dimensions (mm): | 95*80mm  |
| Net weight:      | 259g     |
| Maximum power:   | 26W      |
| Working Voltage: | 7.4V     |

Before using the body slimming and firming device, please read the following information carefully regarding basic safety rules. When using the device, please follow the following precautions.

#### Device description:

1. Charging port
2. RF
3. Cavitation
4. Ring of red light
5. Function change button
6. On/Off button
7. Control panel

#### Precautions:

1. If your skin is very sensitive, shorten the action time for the first use.
2. Avoid places where the device may fall into the sink or bathtub. If the device has been dropped into water, do not reach for it before unplugging the power cord.
3. Do not use any charger other than the original one to charge the device.
4. Do not use the device if it is damaged, wet or not working properly.
5. Do not use the device if you are pregnant or menstruating, have a pacemaker, suffer from epilepsy, heart problems or severe skin diseases.
6. Do not use the device on an area of skin that has undergone plastic surgery.
7. Do not use the device on the mouth, eyeballs or throat.
8. The device is intended for use on healthy skin. Do not use on healing or open wounds.
9. Stop using the device if you feel uncomfortable and consult a doctor immediately.
10. Never attempt to disassemble the device.
11. Keep them out of the reach of children.
12. Keep the device away from fire.
13. Keep the device away from moisture. Do not immerse it in water.
14. Do not drop or hit it.
15. Unplug the power adapter after charging the device.

#### Device features:

1. Modeling the silhouette lines
2. Reduction of excess fat tissue
3. Eliminating cellulite
4. Firming, increasing elasticity and improving the smoothness of the body
5. Skin lifting

#### Radio waves (RF):

Radio frequency (RF) is a proven technology used in the medical and cosmetic industries. At a strictly controlled temperature, the magnetic pulse acts directly on the cells, generating thermal energy. This makes body rejuvenation using radio waves (RF) the best non-surgical facelift for your skin. The use of radio waves will help you effectively shape your figure. You will even eliminate stubborn excess fat tissue and cellulite. Radio wave stimulation reaches both superficial and deep skin tissues. This will ensure that your skin is properly oxygenated, nourished and has proper microcirculation. As a result, while reducing unnecessary fat tissue, you also maintain body firmness at all times.

#### Cavitation (CAVI):

Designed for all skin types. Thanks to a strictly controlled frequency, cavitation penetrates the skin and vibrates at a very high frequency (800KHz). Functions: mechanical, chemical and heating massage even deep skin cells, effectively shrink fat cells and increase the permeability of the cell membrane. This way you will improve circulation, cleanse the lymphatic system and repair dull skin tone. Thanks to cavitation, you will achieve delightful skin cleansing effects. They are visible much faster than with other treatments, and the result of a beautiful, smooth and youthful body lasts for a long time. At the same time, cavitation will help your skin increase the absorption of body care products. At the same time, you can be sure that the cavitation treatment is painless and completely safe.

**EMS:**

EMS electrostimulation stimulates the muscles to work. When using the device, muscle contraction is induced, identical to that experienced during physical exercise. However, 50 times more energy is used than during regular exercise. Such muscle stimulation will allow you to maximally engage your muscles and increase the level of protein and elastic fibers. As a result, you will radically reduce the level of excess fat tissue and cellulite. EMS will restore elasticity and elasticity to your skin and ensure body firming.

**Red photon light:**

Significantly stimulates the production of collagen proteins. Massage using red photon light will allow you to visibly rejuvenate your body by improving its smoothness and elasticity. While reducing fat tissue, you will radically stimulate your blood circulation and thus avoid skin sagging. Using red photon light will give your body vitality, firmness and health. This type of light is perfect for reducing scars, stretch marks and cellulite, as well as reducing and alleviating inflammation.

**Landing:**

1. Before first use, charge the device for at least 3 hours to activate the internal lithium battery.
2. No power: When the battery is low, the device will beep 10 times, then turns off.
3. When the device is placed on the charging cradle, the charging indicator will appear on the screen and the battery icon will gradually change over and over from empty to full. The backlight turns off after 30 seconds of charging.
4. When charging is complete, the charging indicator disappears and the battery icon displays as full.

**Before use:**

Confirmation that the device is ready for cavitation treatment: Lay the cavitation head flat (head up), select the CAVI function and drop 1-2 drops of water on the head. The water droplets will vibrate and be atomized by cavitation.

**Allergy test:**

1. Clean one area of your body on your arm.
2. Apply a skin care product (gel/lotion) to cleansed skin, turn on the device and select your preferred function. Massage this area for 2-3 minutes according to the instructions in this manual and then leave this area of your body without any additional treatments.
3. Wait and check for any irritation or redness within 24 hours.

**Function combination operation manual:**

1. Cleanse your skin.
2. Apply a care product (gel/lotion) evenly to the skin.
3. Turn on the device.
4. Press MODE once to start the COMBI function. Select your preferred intensity level, RF, EMS cavitation will run for 10 seconds and change over and over.
5. Tip: For best results, always massage your body in small circles. If necessary, add gel/lotion during treatment.

**Radio Frequency (RF) Operation Manual**

1. Cleanse your skin.
2. Apply skin care product (gel/lotion) evenly to the skin.
3. Turn on the device.
4. Press MODE twice to start the RF function. Select your preferred intensity level.
5. Tip: For best results, always massage your body in small circles. If necessary, add gel/lotion during treatment.

**Cavitation Instruction Manual (CAVI)**

1. Cleanse your skin.
2. Apply a care product (gel/lotion) evenly to the skin.
3. Turn on the device.
4. Press MODE three times to start the CAVI function. Select your preferred intensity level.
5. Tip: For best results, always massage your body in small circles. If necessary, add gel/lotion during treatment.

**EMS User Manual**

1. Cleanse your skin.
2. Apply the care product (gel/lotion) evenly to the skin.
3. Turn on the device.
4. Press MODE four times to start the EMS function. Select your preferred intensity level.
5. Tip: For best results, always massage your body in small circles. If necessary, add gel/lotion during treatment.

#### Massage instructions:

##### Both arms

Massage both arms evenly, spending the same amount of time on each of them. Massage your arms starting from the elbows towards the armpits. Massage your shoulders in small circles. This form of massage will relieve the tension of your muscles, help you shape your body and maintain its elasticity.

##### Waist and belly

Massage your waist up and down. Then massage your waist from outside to inside. Perform a circular massage around the navel (clockwise). Then massage your belly in large circles (also clockwise). A waist and abdominal massage performed in this way will help you quickly break down even the most stubborn fat tissue, making it easy to get rid of its excess.

##### Leg massage

Massage from the ankles up. For a moment, massage your legs in straight lines upwards (always towards the heart). Then massage in small circles - from the bottom of the leg up (towards the heart). By focusing on body parts particularly susceptible to cellulite, you will effectively get rid of fat tissue. At the same time, you can be sure that your body will not lose its firmness.

##### Hips and buttocks massage

Massage both parts of the body with movements from bottom to top. This massage is used to lift the buttocks and firm the skin.

##### Foot massage

Massage the bottom of your feet, spending the same amount of time on each of them. Massage your feet in a spiral motion. This way you will be free from foot swelling. You will improve circulation, feel relaxed and light - you will get the effect of rested legs.

#### Warnings:

1. Treatments should not exceed 10 minutes for one body part and should be performed a maximum of 3 times a week. The time between each treatment should include at least one day without treatment.
2. When using the COMBI, RF, CAVI or EMS function, always maintain an appropriate amount of care product (lotion/gel) between the massage heads and the skin, especially when using the RF function.
3. Always move the massage heads slowly during the treatment. Leaving the device in one position for a long time may cause skin irritation.

#### After using the device:

1. Turn off your device.
2. Wipe off any remaining skin care product used during the massage on the massage heads with a damp tissue or cloth. (Do not clean it in water and/or any chemical cleaning agent.)
3. Wipe and clean the device with a dry tissue or cloth. Store them in a clean and dry place.

More information about the company and products [page 2 - f2]:

- [www.beautify.eu](http://www.beautify.eu)
- [www.youtube.com/beautifyeu](http://www.youtube.com/beautifyeu)

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**[EN] Rules for handling electrical and electronic equipment:** Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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**[EN] Summary declaration of conformity:** The manufacturer hereby declares under his sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of legislation. Member States in the field of electromagnetic compatibility and Directive 2011/65/EC of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment and Commission Delegated Directive (EU) 2015/863 of 2011. On 31 March 2015, amends Annex II to Directive 2011/65/EU of the European Parliament and of the Council on the list of restricted substances.  
Model: Booster Body