

BEAUTIFLY

Blink
Blush
series

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Eye and lips massager

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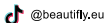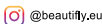f1

f2

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- www.beautifly.eu
- www.youtube.com/beautiflyeu

Dimensions:	119x23x23mm
Net weight:	32g
Temperature:	40OC±5OC
Material:	ABS+PC
Power:	≤1W
Vibration frequency:	12000rpm
Working current:	20-500mA
Working Voltage:	3.7V
Charging voltage:	5V
Charging time:	2h
Battery use time:	2h

Security:

Before using the product, read this user manual and follow the recommendations contained therein when using the device.

Device description [f1]:

1. Head with red-blue LED light (level 1 is red light, level 2 is blue light)
2. Main body
3. USB charging port
4. Massage head
5. On/Off button to switch mode
6. Sensor board
7. Plug

A description of the action:

Hold the On/Off button for 2s to turn the device on or off, short press again switches the red and blue photon light. The red LED will flash when the device is turned on while the massage head is heating up. After reaching the full temperature of 40oC, the diode will glow steadily. The Blink device will vibrate when the massage head and sensing area touch the human body simultaneously.

Area of use:

1. Eye area: Combine with eye care cream/gel, moving evenly and gently massaging from side to side to lift, smooth and deeply nourish the eye area.
2. Forehead (lion's wrinkle): massage upwards from the eyebrow to the hairline, high-frequency micro-vibration massage, acts on skin cells, promotes the penetration of nutrient ions, relieves the lion's frown and relaxes the grip.
3. Cheeks: massaging the cheeks upwards and outwards, accelerates blood circulation, stimulates skin cells, breaks down excess fat on the face, removes fine wrinkles, firms and slims the face.
4. Mouth area: to lift the corners of the lips, improve their blood supply and smooth them, massage from the bottom up.
5. Neck skin: massage from the bottom up. The skin will completely absorb the essence and remain deeply moisturized.

Function description:

- **Intellisense** function, i.e. touch control. If it is not used for 2 minutes, it can be turned off automatically.
- heating function, high-frequency micro-vibration, photorejuvenation and ionization functions.
- high-frequency sonic massage 12,000 times/m, stimulates circulation, accelerates collagen production, restores elastic fibers.
- a comfortable massage head for thermal care maintains a temperature of $40^{\circ}\text{C} \pm 5^{\circ}\text{C}$, opens pores and refreshes the skin, keeping the skin firm, radiant and elastic.
- blue light and red photon light increases skin permeability. Active ingredients penetrate deeper parts of the skin, which can effectively increase its elasticity and reduce melanin around the eyes.
- elliptical arc design of the massage head, closely fits the skin, lift and massage the skin all around, including the corners of the nose and mouth.

Instructions for use:

1. Apply the care product, skin care essence to your face.
2. According to your skin care needs, short press the On/Off button to switch red and blue light (default as red light when you turn on the device).
3. Hold the sensor plate with one hand and touch the skin with the massage head.
4. After use, clean the massage head and replace the cap.

Precautions:

Test the device on your wrist before use. Stop using the massager if it does not produce heat, so as not to damage the skin. If the power cord is damaged, do not use it. The device should not be used in people with: infectious diseases, heart diseases, high fever, malignant tumor, acute inflammation, heart disease, skin diseases, pregnant/breastfeeding. The device uses USB charging (5V), it needs 2 hours to fully charge. When using the device, do not leave the massage head in the same place for too long. Do not disassemble the device yourself. If any fault or improper functioning occurs, discontinue use and return to the designated service center. The device is not waterproof, do not immerse it in water or clean it with corrosive substances. Before cleaning, turn off the power and wait until the device cools down. Clean the head using a cotton towel or tissue. The device should be stored in a dry and ventilated place, avoid direct sunlight. The device is not a toy. Please keep out of the reach of children.

More information about the company and products [f2]:

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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